


Acamprosate

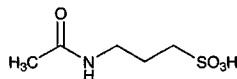

Molecular formula: C₅H₁₁NO₄S

Molecular weight: 181.21

CAS Registry No.: 77337-76-9

SAMPLE

Matrix: blood, urine

Sample preparation: Add 1 mL whole blood or urine to Toxi-Tube A (Toxi-Lab, Irvine CA), add 3 mL water, mix by gentle inversion for 5 min, centrifuge at 1500 g for 5 min. Remove the organic layer and evaporate it to dryness under a stream of nitrogen at 40°, reconstitute the residue with 50 µL MeCN:water 50:50, vortex for 10 s, centrifuge at 7500 g for 2 min, inject a 10 (urine) or 30 (blood) µL aliquot. (The detector wavelength shown is the wavelength of maximum absorbance. This will not necessarily be the optimal wavelength for the separation. Multiple wavelengths from 200-350 nm can be scanned using a diode-array detector. Otherwise, 220 nm may be a reasonable choice for initial work. Matrix may interfere.)

HPLC VARIABLES

Guard column: 20 mm long Symmetry C18

Column: 250 × 4.6 5 µm Symmetry C8 (Waters)

Mobile phase: Gradient. A was 50 mM pH 3.8 sodium phosphate buffer. B was MeCN. A: B 85:15 for 6.5 min, 65:35 for 18.5 min, 20:80 for 3 min (step gradient), re-equilibrate at initial conditions for 7 min.

Column temperature: 30

Flow rate: 1 for 6.5 min, to 1.5 over 18.5 min, maintain at 1.5 for 3 min (re-equilibrate at 1.5 mL/min)

Injection volume: 10-30

Detector: UV 200.5

CHROMATOGRAM

Retention time: 2.97

KEY WORDS

whole blood

REFERENCE

Gaillard, Y.; Pépin, G. Use of high-performance liquid chromatography with photodiode-array UV detection for the creation of a 600-compound library. Application to forensic toxicology, *J. Chromatogr. A*, **1997**, 763, 149-163.

Acebutolol


Molecular formula: C₁₈H₂₈N₂O₄

Molecular weight: 336.43

CAS Registry No.: 37517-30-9

Merck Index: 16

Lednicer No.: 2 109

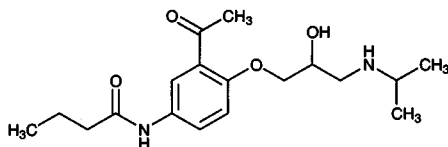

SAMPLE

Matrix: blood

Sample preparation: 200 μ L Serum + 50 μ L water + 500 μ L MeCN, vortex for 1 min. Centrifuge at 6000 rpm for 10 min. Evaporate the supernatant to 200 μ L at 40° under a stream of nitrogen, vortex for 30 s. Inject a 30-80 μ L aliquot.

HPLC VARIABLES

Column: 100 \times 8.0 4 μ m Radial-pak Novapak C18

Mobile phase: MeCN:buffer 14:86 (Buffer was 2 g citric acid, 2 g sodium acetate, and 1 mL triethylamine in 1 L water.)

Flow rate: 2.5

Injection volume: 30-80

Detector: F ex 330 em 440

CHROMATOGRAM

Retention time: 7.4

OTHER SUBSTANCES

Extracted: norfloxacin, pefloxacin

Simultaneous: ciprofloxacin, lomefloxacin, ofloxacin

KEY WORDS

acebutolol is IS; serum

REFERENCE

Abanmi,N.; Zaghlood,I.; El Sayed,N.; al-Khamis,K.I. Determination of pefloxacin and its main active metabolite in human serum by high-performance liquid chromatography, *Ther.Drug Monit.*, **1996**, *18*, 158-163.

SAMPLE

Matrix: blood

Sample preparation: Add 1 mL 67 mM pH 7.4 phosphate buffer, 200 μ L 1 M NaOH, and 6 mL MTBE to 1 mL plasma. Shake for 10 min, centrifuge at 1300 g at 4° for 5 min, freeze the aqueous layer in acetone/dry ice. Add 200 μ L 10 mM HCl to the organic layer, shake for 10 min, centrifuge at 1300 g for 5 min. Freeze the aqueous layer, discard the organic phase, eliminate traces of the organic layer using a stream of cold air over the aqueous layer for 3-4 min. Thaw the aqueous layer, inject a 20 μ L aliquot.

HPLC VARIABLES

Column: 150 \times 4.6 5 μ m Spherisorb hexyl

Mobile phase: MeCN:15 mM pH 3.5 KH₂PO₄ containing 0.05% triethylamine 45:55

Flow rate: 1

Injection volume: 20

Detector: UV 238

CHROMATOGRAM

Retention time: 3.3

Internal standard: acebutolol

OTHER SUBSTANCES

Extracted: celiprolol

KEY WORDS

plasma; acebutolol is IS

REFERENCE

Verbesselt,R.; Zugravu,A.; Tjandramaga,T.B.; De Schepper,P.J. Liquid chromatographic determination of total celiprolol or (S)-celiprolol and (R)-celiprolol simultaneously in human plasma, *J.Chromatogr.B*, **1996**, 683, 231-236.

SAMPLE

Matrix: blood

Sample preparation: 500 μ L Plasma + 200 μ L 1 M NaOH + 5 mL chloroform, vortex for 30 s, centrifuge at 1800 g for 5 min. Remove the organic layer and evaporate it to dryness under reduced pressure, reconstitute the residue in 200 μ L 0.05% S-(+)-1-(1-naphthyl)ethylisocyanate in chloroform, vortex for 30 s, evaporate to dryness under reduced pressure, reconstitute with 200 μ L chloroform, inject a 50-175 μ L aliquot.

HPLC VARIABLES

Column: 250 \times 4.6 Partisil 5 silica

Mobile phase: Hexane:chloroform:MeOH 60:38:2

Flow rate: 2

Injection volume: 50-175

Detector: F ex 220 em 345

CHROMATOGRAM

Retention time: 14.9 (R), 16.5 (S)

Internal standard: (\pm)-acebutolol

OTHER SUBSTANCES

Extracted: tocinide

KEY WORDS

plasma; derivatization; chiral; acebutolol is IS

REFERENCE

Carr,R.A.; Foster,R.T.; Freitag,D.; Pasutto,F.M. Stereospecific high-performance liquid chromatographic determination of tocinide, *J.Chromatogr.*, **1991**, 566, 155-162.

SAMPLE

Matrix: blood

Sample preparation: 500 μ L Plasma + 100 μ L 70 mM pH 7 phosphate buffer + 4 mL chloroform:isopentyl alcohol:diethyl ether 71.25:3.75:25, vortex for 30 s, centrifuge at 1800 g for 5 min. Remove the organic layer and evaporate it to dryness under reduced pressure, reconstitute the residue in 100 μ L chloroform:triethylamine 100:1, add 100 μ L 1% (S)-(+)-1-(1-naphthyl)ethyl isocyanate in chloroform, after 1 min add 50 μ L 2% ethylchloroformate in chloroform, after 30 s add 50 μ L 2.5% ethanolamine in chloroform, inject a 20-125 μ L aliquot.

HPLC VARIABLES

Column: 100 \times 8 4 μ m Nova-Pak silica Radial Pak

Mobile phase: Hexane:chloroform:MeOH 64.5:33:2.5

Flow rate: 2

Injection volume: 20-125

Detector: F ex 245 em 420

CHROMATOGRAM

Retention time: 8.5, 9.5 (enantiomers)

Internal standard: acebutolol

OTHER SUBSTANCES

Extracted: lomefloxacin (F ex 280 em 470)

KEY WORDS

plasma; derivatization; chiral; normal phase; acebutolol is IS

REFERENCE

Foster,R.T.; Carr,R.A.; Pasutto,F.M.; Longstreth,J.A. Stereospecific high-performance liquid chromatographic assay of lomefloxacin in human plasma, *J.Pharm.Biomed.Anal.*, **1995**, 13, 1243–1248.

SAMPLE

Matrix: blood

Sample preparation: 200 μ L Plasma + 50 μ L 50 mM pH 7.4 phosphate buffer + 500 μ L 2% zinc sulfate in MeOH:water 50:50, mix, centrifuge at 13000 rpm for 5 min, inject an aliquot.

HPLC VARIABLES

Guard column: 40 \times 4.6 SynChropak bulk support (Knauer)

Column: 120 \times 4.6 5 μ m Spherisorb ODS1 C18

Mobile phase: MeCN:MeOH:pH 4.5 acetate buffer (ratio not given)

Flow rate: 1

Detector: UV 233

CHROMATOGRAM

Retention time: 4.68

OTHER SUBSTANCES

Extracted: cyclopropane carboxylic acid ester prodrug

KEY WORDS

plasma

REFERENCE

Hovgaard,L.; Brondsted,H.; Buur,A.; Bundgaard,H. Drug delivery studies in Caco-2 monolayers. Synthesis, hydrolysis, and transport of O-cyclopropane carboxylic acid ester prodrugs of various β -blocking agents, *Pharm.Res.*, **1995**, 12, 387–392.

SAMPLE

Matrix: blood

Sample preparation: 2 mL Whole blood or plasma + 2 mL buffer + 5 mL chloroform: isopropanol:n-heptane 60:14:26, shake gently horizontally for 10 min, centrifuge at 2800 g for 10 min. Remove the lower organic layer and evaporate it to dryness under vacuum at 45°, reconstitute the residue in 100 μ L mobile phase, centrifuge at 2800 g for 5 min, inject a 50 μ L aliquot of the supernatant. (Buffer was saturated ammonium chloride solution 25% diluted with water, adjusted to pH 9.5 with 25% ammonia solution.)

HPLC VARIABLES

Column: 300 \times 3.9 4 μ m NovaPack C18

Mobile phase: MeOH:THF:buffer 65:5:30 (Buffer was 0.68 g/L (10 mM (sic)) KH_2PO_4 adjusted to pH 2.6 with concentrated orthophosphoric acid.) (At the end of each session wash the column with water for 1 h and MeOH for 1 h, re-equilibrate for 30 min.)

Column temperature: 30

Flow rate: 0.8

Injection volume: 50

Detector: UV 235

CHROMATOGRAM

Retention time: 3.80

Limit of detection: <120 ng/mL

KEY WORDS

whole blood; plasma; interferences may occur—compounds(all of which are extracted) elute in this order tenoxicam; iproniazid; methocarbamol; methotrexate; caffeine; nialamide; colchicine; cytarabine; benzoylecgonine; acetaminophen; diazoxide; dacarbazine; sulfinpyrazole; flumazenil; sulpride; morphine; atenolol; toloxatone; terbutaline; albuterol; phenobarbital; ranitidine; tiapride; phenol; chlormezanone; aspirin; metformin; ritodrine; codeine; sultopride; amisulpride; naltrexone; lisinopril; benzocaine; nizatidine; nalorphine; mephenesin; naloxone; sotalol; carteolol; procainamide; carbamazepine; bromazepam; nalbuphine; nadolol; procarbazine; dihydralazine; omeprazole; strychnine; acebutolol; glutethimide; chlorpropamide; glipizide; triazolam; prazosin; flunitrazepam; clonazepam; metoclopramide; melfalan; estazolam; tolbutamide; ephedrine; clonidine; pindolol; clobazam; minoxidil; disopyramide; nitrazepam; dextromethorphan; tofisopam; zopiclone; debrisoquine; sulindac; alprazolam; cycloguanil; lorazepam; methaqualone; ketamine; piroxicam; metoprolol; nifedipine; quinine; mephentermine; prilocaine; pentazocine; oxazepam; tiaprofenic acid; quinidine; celiprolol; ajmaline; yohimbine; lidocaine; secobarbital; viloxazine; mepivacaine; meperidine; doxylamine; labetalol; temazepam; amodiaquine; benperidol; droperidol; hydroxychloroquine; zolpidem; ketoprofen; alminoprofen; cicletanine; moclobemide; chloroquine; cocaine; timolol; nomifensine; ticlopidine; acenocoumarol; videsine; mexiletine; dipyridamole; trazodone; pipamperone; pyrimethamine; benazepril; vincristine; metapramine; chlordiazepoxide; oxprenolol; warfarin; clorazepate; flecainide; phenacyclidine; thiopental; fenfluramine; metipranolol; triprolidine; naproxen; buprenorphine; verapamil; buspirone; tianeptine; midazolam; bupivacaine; carbinoxamine; loprazolam; cetirizine; chlorpheniramine; moperone; cibenzoline; medifoxamine; astemizole; vinblastine; nicardipine; bisoprolol; diltiazem; glibornuride; reserpine; aconitine; nitrendipine; diazepam; mianserin; ramipril; haloperidol; tetracaine; alprenolol; aceprometazine; glibenclamide; chlorophenacinone; doxepin; nimodipine; diphenhydramine; cyclizine; histapyrodine; phenylbutazone; demexiptiline; clozapine; proguanil; trifluoperidol; medazepam; cyamemazine; bumadizone; suriclone; propranolol; acepromazine; dothiepin; dextromoramide; fenoprofen; dextropropoxyphene; loxapine; betaxolol; propafenone; promethazine; thioproperazine; methadone; amoxapine; quinupramine; opipramol; cyproheptadine; brompheniramine; mefenidramine; protriptyline; flurbiprofen; tetrazepam; zorubicin; prazepam; alimemazine; loperamide; imipramine; desipramine; levomepromazine; hydroxyzine; niflumic acid; penbutolol; fluvoxamine; pimozide; daunorubicin; indomethacin; maprotiline; tropatenine; etodolac; fluoxetine; amitriptyline; nortriptyline; tiocloamarol; diclofenac; mefloquine; trimipramine; chlorambucil; lidoflazine; ibuprofen; floctafenine; alpidem; loratadine; chlorpromazine; clomipramine; carpipramine; thioridazine; fentiazac; clemastine; mefenamic acid; fluphenazine; prochlorperazine; penfluridol; bepridil; terfenadine; trifluoperazine

REFERENCE

Tracqui,A.; Kintz,P.; Mangin,P. Systematic toxicological analysis using HPLC/DAD, *J.Forensic Sci.*, **1995**, *40*, 254–262.

SAMPLE

Matrix: blood, urine

Sample preparation: Plasma. 1 mL Plasma + 100 μ L 50 μ g/mL pindolol in MeOH + 150 μ L 1 M NaOH + 5 mL diethyl ether, vortex for 30 s, centrifuge at 1800 g for 5 min. Remove the organic layer and evaporate it to dryness under vacuum, add 200 μ L 0.1%

S-(+)-1-(1-naphthyl)ethylisocyanate in chloroform, mix for 30 s, inject a 15-200 μL aliquot. Urine. Dilute 100 fold with water, proceed as above.

HPLC VARIABLES

Column: 250 mm long 5 μm Partisil silica
Mobile phase: Hexane:chloroform:MeOH 63:35:2
Flow rate: 2
Injection volume: 15-200
Detector: F ex 220 em 389

CHROMATOGRAM

Retention time: 12 (R), 13 (S)
Internal standard: pindolol (6,7)
Limit of detection: 1 ng/mL

OTHER SUBSTANCES

Also analyzed: atenolol, nadolol, propranolol, sotalol, toliprolol, tocainide

KEY WORDS

plasma; chiral; derivatization; normal phase

REFERENCE

Piquette-Miller, M.; Foster, R.T.; Pasutto, F.M.; Jamali, F. Stereospecific high-performance liquid chromatographic assay of acebutolol in human plasma and urine, *J.Chromatogr.*, **1990**, 526, 129-137.

SAMPLE

Matrix: blood, urine

Sample preparation: Serum. Condition a 3 mL Supelclean LC-18 SPE cartridge (Supelco) with MeOH and water. Hydrolyze 900 μL serum with β -glucuronidase (EC 3.2.1.31 type H-1 from *Helix pomatia*) at 60° for 1 h, add 500 μL (?) MeOH, centrifuge at 2000 g, add the supernatant to the SPE cartridge, wash with 1 mL water, dry under vacuum, elute with 2 mL MeOH:water 90:10, filter, inject an aliquot. Urine. 900 μL Urine + 500 μL MeOH, filter, inject an aliquot of the filtrate.

HPLC VARIABLES

Guard column: 10 \times 4.6 5 μm HP C18
Column: 150 \times 4.6 5 μm C8P-50 (Asahipak)
Mobile phase: Gradient. MeOH:buffer 30:70 for 4 min, to 45:55 over 6 min, to 50:50 over 2 min, to 60:40 over 2 min, re-equilibrate at initial conditions for 10 min. (Prepare buffer by mixing 100 mM NaH_2PO_4 and 100 mM Na_2HPO_4 to achieve a pH of 7.0 and adding 10 mM N-cetyl-N,N,N-trimethylammonium bromide.)
Injection volume: 20
Detector: UV 260

CHROMATOGRAM

Retention time: 5

OTHER SUBSTANCES

Extracted: alprenolol, atenolol, metoprolol, oxprenolol, propranolol

KEY WORDS

serum; comparison with CE; SPE

REFERENCE

Lukkari, P.; Sirén, H. Ion-pair chromatography and micellar electrokinetic capillary chromatography in analyzing β -adrenergic blocking agents from human biological fluids, *J.Chromatogr.A*, **1995**, 717, 211-217.

SAMPLE

Matrix: blood, urine

Sample preparation: Add 1 mL whole blood or urine to Toxi-Tube A (Toxi-Lab, Irvine CA), add 3 mL water, mix by gentle inversion for 5 min, centrifuge at 1500 g for 5 min. Remove the organic layer and evaporate it to dryness under a stream of nitrogen at 40°, reconstitute the residue with 50 μ L MeCN:water 50:50, vortex for 10 s, centrifuge at 7500 g for 2 min, inject a 10 (urine) or 30 (blood) μ L aliquot. (The detector wavelength shown is the wavelength of maximum absorbance. This will not necessarily be the optimal wavelength for the separation. Multiple wavelengths from 200-350 nm can be scanned using a diode-array detector. Otherwise, 220 nm may be a reasonable choice for initial work. Matrix may interfere.)

HPLC VARIABLES

Guard column: 20 mm long Symmetry C18

Column: 250 \times 4.6 5 μ m Symmetry C8 (Waters)

Mobile phase: Gradient. A was 50 mM pH 3.8 sodium phosphate buffer. B was MeCN. A: B 85:15 for 6.5 min, 65:35 for 18.5 min, 20:80 for 3 min (step gradient), re-equilibrate at initial conditions for 7 min.

Column temperature: 30

Flow rate: 1 for 6.5 min, to 1.5 over 18.5 min, maintain at 1.5 for 3 min (re-equilibrate at 1.5 mL/min)

Injection volume: 10-30

Detector: UV 233.4

CHROMATOGRAM

Retention time: 10.233

KEY WORDS

whole blood

REFERENCE

Gaillard, Y.; Pépin, G. Use of high-performance liquid chromatography with photodiode-array UV detection for the creation of a 600-compound library. Application to forensic toxicology, *J. Chromatogr. A*, **1997**, 763, 149-163.

SAMPLE

Matrix: bulk

Sample preparation: Dissolve 10 μ mole compound (as free base or hydrochloride) in 500 μ L MeCN, add 250 μ L 5% sodium carbonate (for hydrochlorides only), add 500 μ L 100 mM reagent in MeCN, vortex for 1 min, heat at 60° for 2 h, add 100 μ mole L-proline, heat at 60° for 30 min. Remove a 100 μ L aliquot and dilute it with mobile phase, neutralize with acetic acid, inject a 10 μ L aliquot. Prepare the reagent ((R,R)-N-(3,5-dinitrobenzoyl)-2-aminocyclohexylisothiocyanate) as follows. Add 0.7 mL carbon disulfide to 6 mL (1R,2R)-(-)-1,2-diaminocyclohexane, 12 mL water, and 12 mL EtOH, heat the oil bath to 80°, add 2.8 mL carbon disulfide dropwise (making sure that the product does not start to precipitate), when addition is complete reflux for 1 h, acidify with 500 μ L 5 M HCl, reflux for 12 h, cool, filter, wash the solid with a little cold EtOH to give trans-4,5-tetramethyleneimidazolidine-2-thione as a white fluffy solid (mp 148-150°) (Tetrahedron 1993, 49, 4419). Stir 7.97 g 3,5-dinitrobenzoyl chloride in 30 mL dichloroethane at 50°, add a solution of 6 g trans-4,5-tetramethyleneimidazolidine-2-thione in 120 mL dichloroethane containing a catalytic amount of 4-(dimethylamino)pyridine over 15 min, reflux for 2 h, remove the crystals of (R,R)-N-(3,5-dinitrobenzoyl)-2-aminocyclohexylisothiocyanate by filtration, evaporate the filtrate to dryness and dissolve the residue in 60 mL dichloroethane, reflux for 16 h to obtain more (R,R)-N-(3,5-dinitrobenzoyl)-2-aminocyclohexylisothiocyanate (mp >250°, $[\alpha]_{D46} = -133^\circ$ (c = 1) in MeCN).

HPLC VARIABLES

Column: 125 \times 4 5 μ m Lichrospher 60 RP Select B

Mobile phase: MeCN:20 mM ammonium acetate 55:45

Flow rate: 1

Injection volume: 10

Detector: UV 254

CHROMATOGRAM

Retention time: k' 3.10, k' 4.20 (enantiomers)

OTHER SUBSTANCES

Also analyzed: alprenolol, atenolol, carazolol, carvedilol, formoterol, methamphetamine, metipranolol, metoprolol, nifenanol, nitrilo atenolol, oxprenolol, pindolol, propranolol, xamoterol

KEY WORDS

derivatization; chiral

REFERENCE

Kleidermigg,O.P.; Posch,K.; Lindner,W. Synthesis and application of a new isothiocyanate as a chiral derivatizing agent for the indirect resolution of chiral amino alcohols and amines, *J.Chromatogr.A*, **1996**, 729, 33–42.

SAMPLE

Matrix: formulations

Sample preparation: Weigh 10 tablets, powder finely. Weigh accurately powder containing 10 mg nifedipine, dissolve in MeOH and make up to 50 mL with MeOH. Add 1.6 mg IS, filter through a 45 μ m membrane filter. Inject a 10 μ L aliquot.

HPLC VARIABLES

Column: 250 \times 4.6 10 μ m Lichrosorb RP C18

Mobile phase: MeOH:water 55:45 pH 4.5

Flow rate: 1.0 for 4 min, then 2.0

Injection volume: 10

Detector: UV 260

CHROMATOGRAM

Retention time: 3.31

Internal standard: oxprenolol (4.35)

Limit of detection: 3.01 μ g/mL

OTHER SUBSTANCES

Simultaneous: nifedipine, nifedipine oxidation products

KEY WORDS

comparison with GC and first-derivative spectroscopy; tablets

REFERENCE

el Walily,A.F.M. Analysis of nifedipine--acebutolol hydrochloride binary combination in tablets using UV-derivative spectroscopy, capillary gas chromatography and high performance liquid chromatography, *J.Pharm.Biomed.Anal.*, **1997**, 16, 21–30.

SAMPLE

Matrix: formulations

Sample preparation: Take up in mobile phase, inject an aliquot.

HPLC VARIABLES

Column: 250 \times 4.6 10 μ m LiChrosorb C2

Mobile phase: MeCN:buffer 35:65 (1 mL 100 mM HCl + 1200 mL water + 5.84 g NaCl, mix to dissolve, add 700 mL MeOH, make up to 2 L, apparent pH 4.5.)

Flow rate: 1.2

Injection volume: 20

Detector: UV 254

CHROMATOGRAM

Retention time: 9

OTHER SUBSTANCES

Simultaneous: atenolol, nadolol, alprenolol, metoprolol, oxprenolol, pindolol, practolol, propranolol, sotalol

Interfering: timolol

KEY WORDS

tablets

REFERENCE

Patel,B.R.; Kirschbaum,J.J.; Poet,R.B. High-pressure liquid chromatography of nadolol and other β -adrenergic blocking drugs, *J.Pharm.Sci.*, **1981**, 70, 336-338.

SAMPLE

Matrix: solutions

Sample preparation: Mix 1 mL of an aqueous solution with 1 mL 100 mM nickel sulfate in water, 1 mL 20% aqueous ammonia, and 5 mL chloroform:carbon disulfide 98:2, shake vigorously for 1 min, wash the organic layer with three 2 mL portions of water, filter (phase-separation paper). Evaporate the filtrate to dryness under a stream of nitrogen, reconstitute with 1 mL mobile phase, inject a 10 μ L aliquot. (Copper may also be used with electrochemical detection or UV detection at 270 nm.)

HPLC VARIABLES

Guard column: 30 \times 4 40 μ m LiChrosorb RP-18

Column: 250 \times 4 7 μ m LiChrosorb RP-18

Mobile phase: MeOH:20 mM pH 5.8 sodium acetate buffer 80:20 containing 5 mM lithium perchlorate

Flow rate: 1.5

Injection volume: 10

Detector: UV 325, E, Merck-Clevenot E 230, Model LCC 231 thin-layer electrolytic cell with a glassy carbon electrode at +0.7 V, standard calomel reference electrode

CHROMATOGRAM

Retention time: k' 2.48

Limit of detection: 1 fmole (E), 1 nmole (UV)

OTHER SUBSTANCES

Also analyzed: alprenolol, ephedrine, flecainide, methamphetamine, propranolol

KEY WORDS

derivatization; complexation

REFERENCE

Leroy,P.; Nicolas,A. Determination of secondary amino drugs as their metal dithiocarbamate complexes by reversed-phase high-performance liquid chromatography with electrochemical detection, *J.Chromatogr.*, **1984**, 317, 513-521.

SAMPLE

Matrix: solutions

Sample preparation: Prepare a 10 µg/mL solution in MeOH, inject a 20 µL aliquot.

HPLC VARIABLES

Column: 125 × 4.9 Spherisorb S5W silica

Mobile phase: MeOH containing 10 mM ammonium perchlorate and 1 mL/L 100 mM NaOH in MeOH, pH 6.7

Flow rate: 2

Injection volume: 20

Detector: E, LeCarbone, V25 glassy carbon electrode, + 1.2 V

CHROMATOGRAM

Retention time: 2.2

OTHER SUBSTANCES

Also analyzed: acepromazine, acetophenazine, N-acetylprocainamide, albuterol, alprenolol, amethocaine, amiodarone, amitriptyline, antazoline, atenolol, azacyclonal, bamethan, benactyzine, benperidol, benzethidine, benzocaine, benzocetamine, benzphetamine, benzquinamide, bromhexine, bromodiphenhydramine, bromperidol, brompheniramine, brompromazine, buclizine, bufotenine, bupivacaine, buprenorphine, butacaine, butethamate, chlorcyclizine, chlorpheniramine, chlorphenoxamine, chlorprenaline, chlorpromazine, chlorprothixene, cimetidine, cinchonidine, cinnarizine, clemastine, clomipramine, clonidine, cocaine, cyclazocine, cyclizine, cyclopentamine, cyproheptadine, deserpidine, desipramine, dextromoramide, dextropropoxyphene, dicyclomine, diethylcarbamazine, diethylpropion, diethylthiambutene, dihydroergotamine, dimethindene, dimethothiazine, diphenhydramine, diphenoxylate, dipipanone, diprenorphine, dipyridamole, disopyramide, dothiepin, doxapram, doxepin, doxylamine, droperidol, ephedrine, ergocornine, ergocristine, ergocristinine, ergocryptine, ergometrine, ergosine, ergosinine, ergotamine, ethopropazine, etorphine, etoxeridine, fenethazine, fenfluramine, fenoterol, fentanyl, flavoxate, fluopromazine, flupenthixol, fluphenazine, flurazepam, haloperidol, hydroxyzine, hyoscine, ibogaine, imipramine, indapamine, iprindole, isothipendyl, isoxsuprine, ketanserine, laudanosine, lidocaine, lofepramine, loxapine, maprotiline, mecamlamine, meclophenoxate, meclozine, medazepam, mephentermine, mepivacaine, meptazinol, mepyramine, mesoridazine, metaraminol, methadone, methamphetamine, methapyrilene, methdilazene, methotrimeprazine, methoxamine, methoxyphenamine, methoxypropazine, methylephedrine, methylegonovine, methysergide, metoclopramide, metopimazine, metoprolol, mianserin, morazone, nadolol, nalorphine, naloxone, naphazoline, nicotine, nifedipine, nomifensine, nortriptyline, noscapine, orphenadrine, oxeladin, oxprenolol, oxymetazolin, papaverine, pargyline, pecazine, penbutolol, pentazocine, penthienate, pericyazine, perphenazine, phenadoxone, phenampromide, phenazocine, phenbutrazate, phendimetrazine, phenelzine, phenglutarimide, phenindamine, pheniramine, phenmetrazine, phenomorphan, phenoperidine, phenothiazine, phenoxybenzamine, phentolamine, phenylephrine, phenyltoloxamine, physostigmine, piminodine, pimozide, pindolol, pipamazine, pipazethate, piperacetazine, piperidolate, pipradol, pirenzepine, piritramide, pizotifen, practolol, pramoxine, prazosin, prenylamine, prilocaine, primaquine, proadifen, procainamide, procaine, prochlorperazine, procyclidine, proheptazine, prolintane, promazine, promethazine, pronethalol, properidine, propiomazine, propranolol, prothipendyl, protriptyline, proxymetacaine, pseudoephedrine, pyrimethamine, quinidine, quinine, ranitidine, rescinnamine, sotalol, tacrine, terazosin, terbutaline, terfenadine, thenyldiamine, theophylline, thiethylperazine, thiopropazate, thioproperazine, thioridazine, thiothixene, thonzylamine, timolol, tocanide, tolpropamine, tolycaine, tranlycypromine, trazodone, trifluoperazine, trifluoperidol, trimeperidine, trimeprazine, trimethobenzamide, trimethoprim, trimipramine, tripeleminamine, triprolidine, tryptamine, verapamil, xylometazoline

REFERENCE

Jane, I.; McKinnon, A.; Flanagan, R.J. High-performance liquid chromatographic analysis of basic drugs on silica columns using non-aqueous ionic eluents. II. Application of UV, fluorescence and electrochemical oxidation detection, *J.Chromatogr.*, **1985**, *323*, 191–225.

SAMPLE

Matrix: solutions

Sample preparation: Filter (0.22 μm), inject a 10 μL aliquot.

HPLC VARIABLES

Column: 250 \times 4.6 internal surface reversed-phase silica (Pinkerton) (Regis Chemical)

Mobile phase: Isopropanol:100 mM pH 6.8 KH_2PO_4 10:90

Flow rate: 1

Injection volume: 10

Detector: UV 232-274 (wavelength of maximum absorption used)

CHROMATOGRAM

Retention time: 32.2

OTHER SUBSTANCES

Simultaneous: carteolol, atenolol, metoprolol, oxprenolol, pindolol, alprenolol

REFERENCE

Ohshima,T.; Takagi,K.; Miyamoto,K.-I. High performance liquid chromatographic retention time of β -blockers as an index of pharmacological activity, *J.Liq.Chromatogr.*, **1993**, 16, 3933–3939.

SAMPLE

Matrix: solutions

HPLC VARIABLES

Column: 300 \times 3.9 5 μm Nova-Pak C18

Mobile phase: MeOH:buffer 30:70 (Buffer was pH 4.0 phosphate buffer (ionic strength = 0.1) containing 2.86 mM N,N-dimethyloctylamine, pH readjusted to 4.00 with 85% phosphoric acid.)

Column temperature: 30

Flow rate: 1

Injection volume: 100

Detector: UV 220

CHROMATOGRAM

Retention time: k' 2.54

OTHER SUBSTANCES

Also analyzed: bunitrolol, carazolol, celiprolol, esmolol, mepindolol, metoprolol, timolol

REFERENCE

Hamoir,T.; Verlinden,Y.; Massart,D.L. Reversed-phase liquid chromatography of β -adrenergic blocking drugs in the presence of a tailing suppressor, *J.Chromatogr.Sci.*, **1994**, 32, 14–20.

SAMPLE

Matrix: solutions

HPLC VARIABLES

Column: 62 \times 2 packed with chiral packing (Prepare packing by dissolving 4-chloro-3-methylphenylcarbamate cellulose in THF, coat on Nucleosil 1000-7, dry at 60° for 3 h under reduced pressure.)

Mobile phase: Hexane:isopropanol:diethylamine 90:10:0.1

Flow rate: 0.1

Injection volume: 20

Detector: UV 254

CHROMATOGRAM

Retention time: k' 25.7

KEY WORDS

narrow-bore; chiral; α 1.12

REFERENCE

Chankvetadze,B.; Chankvetadze,L.; Sidamonidze,S.; Yashima,E.; Okamoto,Y. Enantioseparation of some chiral pharmaceuticals using narrow-bore liquid chromatography, *J.Pharm.Biomed.Anal.*, **1995**, 13, 695–699.

SAMPLE

Matrix: solutions

HPLC VARIABLES

Column: 250 × 4.6 Chirex 3022 (Phenomenex)

Mobile phase: Hexane:1,2-dichloroethane:EtOH/trifluoroacetic acid 55:35:10 (EtOH/trifluoroacetic acid was premixed 20:1.)

Flow rate: 0.7-1

Injection volume: 20

Detector: UV 244

KEY WORDS

chiral; α = 1.09 for enantiomers

REFERENCE

Cleveland,T. Pirkle-concept chiral stationary phases for the HPLC separation of pharmaceutical racemates, *J.Liq.Chromatogr.*, **1995**, 18, 649–671.

SAMPLE

Matrix: solutions

Sample preparation: Inject a 20 μ L aliquot of a 1 mg/mL solution.

HPLC VARIABLES

Column: 250 × 4.6 10 μ m Chiralcel OD

Mobile phase: Hexane:EtOH:diethylamine 90:10:0.1

Flow rate: 0.5

Injection volume: 20

Detector: UV 320

CHROMATOGRAM

Retention time: k' 2.77, 3.19 (enantiomers)

KEY WORDS

chiral

REFERENCE

Ekelund,J.; van Arkens,A.; Bronnum-Hansen,K.; Fich,K.; Olsen,L.; Petersen,P.V. Chiral separations of β -blocking drug substances using chiral stationary phases, *J.Chromatogr.A*, **1995**, 708, 253–261.

SAMPLE

Matrix: solutions

HPLC VARIABLES

Column: 150 × 4.6 12 μ m 1-myristoyl-2-[(13-carboxyl)-tridecoyl]-sn-3-glycerophosphocholine chemically bonded to silica (Regis)

Mobile phase: MeCN:100 mM pH 7.0 phosphate buffer 20:80

Flow rate: 1

Detector: UV 254

CHROMATOGRAM

Retention time: k' 4.00

OTHER SUBSTANCES

Also analyzed: alprenolol, antazoline, atenolol, betaxolol, bisoprolol, bopindolol, bupranolol, carteolol, celiprolol, chloropyramine, chlorpheniramine, cicloprolol, cimetidine, cinarizine, cirazoline, clonidine, dilevalol, dimethindene, diphenhydramine, doxazosin, esmolol, famotidine, isothipendyl, ketotifen, metiamide, metoprolol, moxonidine, nadolol, naphazoline, nifenalol, nizatidine, oxprenolol, pheniramine, phentolamine, pindolol, pizotyline (pizotifen), practolol, prazosin, promethazine, propranolol, pyrilamine (mepyramine), ranitidine, roxatidine, sotalol, tiamenidine, timolol, tramazoline, tripeleminamine, triprolidine, tymazoline, UK-14,304

REFERENCE

Kaliszan, R.; Nasal, A.; Turowski, M. Binding site for basic drugs on α_1 -acid glycoprotein as revealed by chemometric analysis of biochromatographic data, *Biomed. Chromatogr.*, **1995**, *9*, 211–215.

SAMPLE

Matrix: solutions

HPLC VARIABLES

Column: 250 × 4.6 5 μ m Supelcosil LC-DP (A) or 250 × 4 5 μ m LiChrospher 100 RP-8 (B)

Mobile phase: MeCN:0.025% phosphoric acid:buffer 25:10:5 (A) or 60:25:15 (B) (Buffer was 9 mL concentrated phosphoric acid and 10 mL triethylamine in 900 mL water, adjust pH to 3.4 with dilute phosphoric acid, make up to 1 L.)

Flow rate: 0.6

Injection volume: 25

Detector: UV 229

CHROMATOGRAM

Retention time: 6.51 (A), 3.72 (B)

OTHER SUBSTANCES

Also analyzed: acepromazine, acetaminophen, acetazolamide, acetophenazine, albuterol, alprazolam, amitriptyline, amobarbital, amoxapine, antipyrine, atenolol, atropine, azatadine, baclofen, benzocaine, bromocriptine, brompheniramine, brotizolam, bupivacaine, buspirone, butabarbital, butalbital, caffeine, carbamazepine, cetirizine, chlorcyclizine, chlordiazepoxide, chlormezanone, chloroquine, chlorpheniramine, chlorpromazine, chlorpropamide, chlorprothixene, chlorthalidone, chlorzoxazone, cimetidine, cisapride, clomipramine, clonazepam, clonidine, clozapine, cocaine, codeine, colchicine, cyclizine, cyclobenzaprine, dantrolene, desipramine, diazepam, diclofenac, diflunisal, diltiazem, diphenhydramine, diphenidol, diphenoxylate, dipyridamole, disopyramide, dobutamine, doxapram, doxepin, droperidol, encainide, ethidium bromide, ethopropazine, fenoprofen, fentanyl, flavoxate, fluoxetine, fluphenazine, flurazepam, flurbiprofen, fluvoxamine, furosemide, glutethimide, glyburide, guaifenesin, haloperidol, homatropine, hydralazine, hydrochlorothiazide, hydrocodone, hydromorphone, hydroxychloroquine, hydroxyzine, ibuprofen, imipramine, indomethacin, ketoconazole, ketoprofen, ketorolac, labetalol, levorphanol, lidocaine, loratadine, lorazepam, lovastatin, loxapine, mazindol, mefenamic acid, meperidine, mephénytoin, mepivacaine, mesoridazine, metaproterenol, methadone, methdilazine, methocarbamol, methotrexate, methotrimeprazine, methoxamine, methyl-dopa, methylphenidate, metoclopramide, metolazone, metoprolol, metronidazole, midazolam, moclobemide, morphine, nadolol, nalbuphine, naloxone, naphazoline, naproxen, nifedipine, nizatidine, norepinephrine, nortriptyline, oxazepam, oxycodone, oxymetazoline, paroxetine, pemoline, pentazocine, pentobarbital, pentoxifylline, perphenazine,

pheniramine, phenobarbital, phenol, phenolphthalein, phentolamine, phenylbutazone, phenyltoloxamine, phenytoin, pimozide, pindolol, piroxicam, pramoxine, prazepam, prazosin, probenecid, procainamide, procaine, prochlorperazine, procyclidine, promazine, promethazine, propafenone, propantheline, propiomazine, propofol, propranolol, protriptyline, quazepam, quinidine, quinine, racemethorphan, ranitidine, remoxipride, risperidone, salicylic acid, scopolamine, secobarbital, sertraline, sotalol, spironolactone, sulfapyrazole, sulindac, temazepam, terbutaline, terfenadine, tetracaine, theophylline, thiethylperazine, thiopental, thioridazine, thiothixene, timolol, tocinide, tolbutamide, tolmetin, trazodone, triamterene, triazolam, trifluoperazine, triflupromazine, trimeprazine, trimethoprim, trimipramine, verapamil, warfarin, xylometazoline, yohimbine, zopiclone

KEY WORDS

also details of plasma extraction

REFERENCE

Koves, E.M. Use of high-performance liquid chromatography-diode array detection in forensic toxicology, *J.Chromatogr.A*, **1995**, 692, 103–119.

SAMPLE

Matrix: solutions

Sample preparation: 100 μ L 55 mM N-Benzoyloxycarbonyl-L-phenylalanine (N-CBZ-L-Phe) in dichloromethane + 100 μ L 14 mM N,N-dimethylaminopyrrole (dimethylaminopyrrole (?)) in dichloromethane + 100 μ L 9-acetylanthracycline in dichloromethane, cool in an ice bath, add 500 μ L of a solution of acebutolol in dichloromethane, add 100 μ L 240 mM dicyclohexylcarbodiimide in dichloromethane, shake mechanically at 0° for 30 min, add 100 μ L 1.06 M acetic anhydride in dichloromethane, shake mechanically at 30° for 15 min, add 1 mL MeOH, mix, inject an aliquot.

HPLC VARIABLES

Guard column: 10 μ m Nova-Pak C18 precolumn

Column: 150 \times 3.9 4 μ m Nova-Pak C18

Mobile phase: MeOH:water 60:40

Flow rate: 1.3

Detector: UV 254

CHROMATOGRAM

Retention time: 14.13 (S), 15.74 (R)

Internal standard: 9-acetylanthracycline (6.52)

KEY WORDS

derivatization; chiral

REFERENCE

Wen, Y.H.; Wu, S.S.; Wu, H.L. Chiral separation of acebutolol by derivatization and high-performance liquid chromatography, *J.Liq.Chromatogr.*, **1995**, 18, 3329–3345.

SAMPLE

Matrix: solutions

Sample preparation: Mix 20 μ L of a 1 mM solution in MeOH or water with 50 μ L pH 8 borate buffer and 50 μ L 18 mM 2-(6-methoxy-2-naphthyl)-1-propyl chloroformate in acetone, vortex, let stand at room temperature for 30 min, add 100 μ L 10 mM trans-4-hydroxy-L-proline in water, mix, let stand for 2 min, add 2 mL dichloromethane, vortex for 30 s. Remove the organic layer and evaporate it to dryness under reduced pressure, reconstitute the residue in 100 μ L mobile phase, inject an aliquot. Prepare 2-(6-methoxy-2-naphthyl)-1-propyl chloroformate as follows. Stir 1.5 mmoles lithium aluminum hydride in THF, slowly add 2 mmoles (S)-naproxen in 20 mL anhydrous THF, reflux for 1 h, evaporate most of the solvent, cautiously add water with stirring, acidify with 6 N HCl,

extract three times with diethyl ether. Combine the organic layers and dry them over anhydrous sodium sulfate, evaporate to dryness, chromatograph on silica gel with dichloromethane:MeOH 100:2 (flash chromatography), evaporate eluate to dryness, dry under vacuum over KOH to give 2-(6-methoxy-2-naphthyl)propanol as a white solid (mp 92-3°). Stir 0.5 mmoles 2-(6-methoxy-2-naphthyl)propanol and 0.5 mmoles triethylamine in 10 mL dry toluene at 0°, add 1 mL 20% phosgene in toluene (Caution! Phosgene is highly toxic, perform reaction in a chemical fume hood!) (Fluka), stir for 4 h, filter, evaporate to dryness under reduced pressure, dry under vacuum to give 2-(6-methoxy-2-naphthyl)-1-propyl chloroformate (mp 60°). Store under vacuum over phosphorus pentoxide at room temperature.)

HPLC VARIABLES

Column: 250 × 4 5 µm Zorbax-SIL

Mobile phase: n-Hexane:isopropanol 100:5

Flow rate: 1.5

Injection volume: 100

Detector: UV 230, F ex 270 em 365

CHROMATOGRAM

Retention time: k' 20.4, 21.9 (enantiomers)

KEY WORDS

derivatization; chiral; normal phase

REFERENCE

Büschges,R.; Linde,H.; Mutschler,E.; Spahn-Langguth,H. Chloroformates and isothiocyanates derived from 2-arylpropionic acids as chiral reagents: synthetic routes and chromatographic behaviour of the derivatives, *J.Chromatogr.A*, **1996**, 725, 323-334.

SAMPLE

Matrix: urine

Sample preparation: Direct injection.

HPLC VARIABLES

Column: 250 × 4 OmniPac PAX-500 (Dionex)

Mobile phase: Gradient. MeCN:100 mM sodium carbonate 4.5:95.5, after 0.1 min MeCN:water 4.5:95.5, inject, stay with this mobile phase for 5 min then go to MeCN:water 67.5:32.5 over 15 min, re-equilibrate for 10 min before next injection.

Flow rate: 1

Detector: UV 254

CHROMATOGRAM

Retention time: 12

KEY WORDS

rat

REFERENCE

Slingsby,R.W.; Rey,M. Determination of pharmaceuticals by multi-phase chromatography: Combined reversed phase and ion exchange in one column, *J.Liq.Chromatogr.*, **1990**, 13, 107-134.

SAMPLE

Matrix: urine

Sample preparation: 1 mL Urine + 10 mg β-glucuronidase/arylsulfatase (Helix pomatia, Sigma), heat at 37° overnight, add an equal volume of buffer, centrifuge at 2000 g for 5 min, inject an aliquot of the supernatant onto column A with mobile phase A and elute to waste. After 2.5 min backflush the contents of column A onto column B with mobile

phase B, monitor the effluent from column B. For gradient elution, after 15 min re-equilibrate both columns for 12.5 min before the next injection. For isocratic elution, remove column A from the circuit after 1.25 min, re-equilibrate column A for 1.5 min. (Buffer was 200 mM boric acid adjusted to pH 9.5 with 5 M NaOH.)

HPLC VARIABLES

Column: A 10 × 4.6 5 μm Spherisorb cyanopropyl; B 250 × 4.6 Capcell Pak C18 UG-120 (Shiseido)

Mobile phase: A water; B Gradient. MeCN:buffer from 3:97 to 30:70 over 30 min, to 40:60 over 8 min (for screening) or isocratic 22:78 (Buffer was 3.4 mL/L phosphoric acid adjusted to pH 3.0 with 5 M NaOH.)

Flow rate: A 1.25; B 1

Injection volume: 100

Detector: UV 235

CHROMATOGRAM

Retention time: 11.5 (gradient), 6 (isocratic)

Limit of detection: 250 ng/mL

OTHER SUBSTANCES

Extracted: (using gradient elution) metabolites, alprenolol, amphetamine, atenolol, bopindolol, codeine, ephedrine, labetalol, metoprolol, morphine, nadolol, oxprenolol, pindolol, propranolol

Interfering: timolol

KEY WORDS

column-switching

REFERENCE

Saarinén, M.T.; Sirén, H.; Riekkola, M.-L. Screening and determination of β-blockers, narcotic analgesics and stimulants in urine by high-performance liquid chromatography with column switching, *J. Chromatogr. B*, **1995**, 664, 341–346.


Aceclofenac

Molecular formula: $C_{16}H_{13}Cl_2NO_4$

Molecular weight: 354.19

CAS Registry No.: 89796-99-6

Merck Index: 19

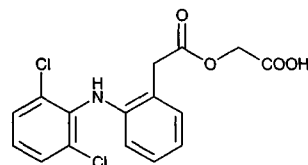

SAMPLE

Matrix: blood

Sample preparation: Acidify plasma, extract with hexane:isopropanol 90:10, evaporate, reconstitute with mobile phase, inject an aliquot.

HPLC VARIABLES

Column: 125×4.6 μ m Spherisorb 5C8

Mobile phase: MeOH:70 mM pH 5.8 phosphate buffer 49:51

Flow rate: 1.7

Detector: UV 282

CHROMATOGRAM

Internal standard: flufenamic acid, niflumic acid

Limit of detection: 3 pmole

OTHER SUBSTANCES

Extracted: metabolites, diclofenac

KEY WORDS

plasma; rat; human; monkey; pharmacokinetics

REFERENCE

Bort,R.; Ponsoda,X.; Carrasco,E.; Gómez-Lechón,M.J.; Castell,J.V. Comparative metabolism of the non-steroidal antiinflammatory drug, aceclofenac, in the rat, monkey, and human, *Drug Metab.Dispos.*, 1996, 24, 969-975.

SAMPLE

Matrix: blood, microsomal incubations, urine

Sample preparation: Deconjugate plasma and urine samples with 50 mU/mL β -glucuronidase and 30 mU/mL arylsulfatase in 100 mM pH 4.5 acetate buffer containing 120 mM NaF, heat at 37° for 4 h. Add an equal volume of MeCN, mix, centrifuge at 9000 rpm for 10 min. Add an equal volume of 100 mM pH 7.4 phosphate buffer to the supernatant, mix, inject an aliquot.

HPLC VARIABLES

Column: 200×4.6 μ m Spherisorb ODS2

Mobile phase: MeCN: 0.02% triethanolamine in 100 mM pH 7.4 phosphate buffer 25:75

Flow rate: 1

Detector: UV 282

CHROMATOGRAM

Retention time: 23

Internal standard: carprofen (11)

OTHER SUBSTANCES

Extracted: metabolites

KEY WORDS

plasma; urine; metabolites

REFERENCE

Bort,R.; Ponsoda,X.; Carrasco,E.; Gómez-Lechón,M.J.; Castell,J.V. Metabolism of aceclofenac in humans, *Drug Metab.Dispos.*, **1996**, 24, 834–841.

SAMPLE

Matrix: cultured hepatocytes, microsomal incubations

Sample preparation: Cultured hepatocytes. Dilute the incubation mixture with an equal volume of MeCN, centrifuge at 9000 rpm for 20 min, dilute 50:50 with 100 mM pH 7.4 phosphate buffer, inject an aliquot. Microsomal incubations. Add an equal volume of cold MeCN to the microsomal incubation, centrifuge, dilute 50:50 with 100 mM pH 7.4 phosphate buffer, inject an aliquot.

HPLC VARIABLES

Column: 200 × 4.6 5 µm Spherisorb ODS2

Mobile phase: MeCN:100 mM phosphate buffer containing 0.02% triethanolamine 25:75

Flow rate: 1

Detector: UV 282

CHROMATOGRAM

Retention time: 23

Internal standard: carprofen (11)

OTHER SUBSTANCES

Extracted: metabolites

KEY WORDS

rat; human; pharmacokinetics

REFERENCE

Bort,R.; Ponsoda,X.; Carrasco,E.; Gómez-Lechón,M.J.; Castell,J.V. Comparative metabolism of the non-steroidal antiinflammatory drug, aceclofenac, in the rat, monkey, and human, *Drug Metab.Dispos.*, **1996**, 24, 969–975.


Acefylline

Molecular formula: C₉H₁₀N₄O₄

Molecular weight: 238.20

CAS Registry No.: 652-37-9, 837-27-4 (sodium salt)

Merck Index: 22

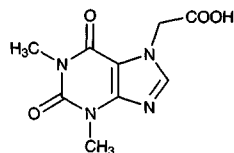

SAMPLE

Matrix: blood, urine

Sample preparation: Add 1 mL whole blood or urine to Toxi-Tube A (Toxi-Lab, Irvine CA), add 3 mL water, mix by gentle inversion for 5 min, centrifuge at 1500 g for 5 min. Remove the organic layer and evaporate it to dryness under a stream of nitrogen at 40°, reconstitute the residue with 50 µL MeCN:water 50:50, vortex for 10 s, centrifuge at 7500 g for 2 min, inject a 10 (urine) or 30 (blood) µL aliquot. (The detector wavelength shown is the wavelength of maximum absorbance. This will not necessarily be the optimal wavelength for the separation. Multiple wavelengths from 200-350 nm can be scanned using a diode-array detector. Otherwise, 220 nm may be a reasonable choice for initial work. Matrix may interfere.)

HPLC VARIABLES

Guard column: 20 mm long Symmetry C18

Column: 250 × 4.6 5 µm Symmetry C8 (Waters)

Mobile phase: Gradient. A was 50 mM pH 3.8 sodium phosphate buffer. B was MeCN. A: B 85:15 for 6.5 min, 65:35 for 18.5 min, 20:80 for 3 min (step gradient), re-equilibrate at initial conditions for 7 min.

Column temperature: 30

Flow rate: 1 for 6.5 min, to 1.5 over 18.5 min, maintain at 1.5 for 3 min (re-equilibrate at 1.5 mL/min)

Injection volume: 10-30

Detector: UV 207.5

CHROMATOGRAM

Retention time: 3.628

KEY WORDS

whole blood

REFERENCE

Gaillard, Y.; Pépin, G. Use of high-performance liquid chromatography with photodiode-array UV detection for the creation of a 600-compound library. Application to forensic toxicology, *J. Chromatogr. A*, **1997**, 763, 149-163.


Acenocoumarol

Molecular formula: C₁₉H₁₅NO₃

Molecular weight: 353.33

CAS Registry No.: 152-72-7

Merck Index: 29

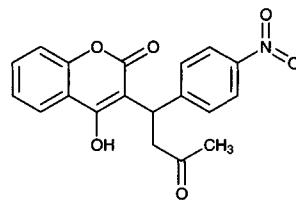

SAMPLE

Matrix: blood

Sample preparation: 2 mL Whole blood or plasma + 2 mL buffer + 5 mL chloroform: isopropanol:n-heptane 60:14:26, shake gently horizontally for 10 min, centrifuge at 2800 g for 10 min. Remove the lower organic layer and evaporate it to dryness under vacuum at 45°, reconstitute the residue in 100 µL mobile phase, centrifuge at 2800 g for 5 min, inject a 50 µL aliquot of the supernatant. (Buffer was saturated ammonium chloride solution 25% diluted with water, adjusted to pH 9.5 with 25% ammonia solution.)

HPLC VARIABLES

Column: 300 × 3.9 4 µm NovaPack C18

Mobile phase: MeOH:THF:buffer 65:5:30 (Buffer was 0.68 g/L (10 mM (sic)) KH₂PO₄ adjusted to pH 2.6 with concentrated orthophosphoric acid.) (At the end of each session wash the column with water for 1 h and MeOH for 1 h, re-equilibrate for 30 min.)

Column temperature: 30

Flow rate: 0.8

Injection volume: 50

Detector: UV 283

CHROMATOGRAM

Retention time: 4.91

Limit of detection: <120 ng/mL

KEY WORDS

whole blood; plasma; interferences may occur—compounds(all of which are extracted) elute in this order tenoxicam; iproniazid; methocarbamol; methotrexate; caffeine; nialamide; colchicine; cytarabine; benzoylecgonine; acetaminophen; diazoxide; dacarbazine; sulfapyrazole; flumazenil; sulpride; morphine; atenolol; tolaxatone; terbutaline; albuterol; phenobarbital; ranitidine; tiapride; phenol; chlormezanone; aspirin; metformin; ritodrine; codeine; sultopride; amisulpride; naltrexone; lisinopril; benzocaine; nizatidine; nalorphine; mephensin; naloxone; sotalol; carteolol; procainamide; carbamazepine; bromazepam; nalbuphine; nadolol; procarbazine; dihydralazine; omeprazole; strychnine; acebutolol; glutethimide; chlorpropamide; glipizide; triazolam; prazosin; flunitrazepam; clonazepam; metoclopramide; melphalan; estazolam; tolbutamide; ephedrine; clonidine; pindolol; clobazam; minoxidil; disopyramide; nitrazepam; dextromethorphan; tofisopam; zopiclone; debrisoquine; sulindac; alprazolam; cycloguanil; lorazepam; methaqualone; ketamine; piroxicam; metoprolol; nifedipine; quinine; mephentermine; prilocaine; pentazocine; oxazepam; tiaprofenic acid; quinidine; celiprolol; ajmaline; yohimbine; lidocaine; secobarbital; viloxazine; mepivacaine; meperidine; doxylamine; labetalol; temazepam; amodiaquine; benperidol; droperidol; hydroxychloroquine; zolpidem; ketoprofen; alminoprofen; cicletanine; moclobemide; chloroquine; cocaine; timolol; nomifensine; ticlopidine; acenocoumarol; vindesine; mexiletine; dipyridamole; trazodone; pipamperone; pyrimethamine; benazepril; vincristine; metapramine; chlordiazepoxide; oxprenolol; warfarin; clorazepate; flecainide; phenacyclidine; thiopental; fenfluramine; metipranolol; triprolidine; naproxen; buprenorphine; verapamil; buspirone; tianeptine; midazolam; bupivacaine; carbinoxamine; loperazolam; cetirizine; chlorpheniramine; moperone; cibenzoline; medifoxamine; astemizole; vinblastine; nicardipine; bisoprolol; diltiazem; glibornuride; reserpine; aconitine; nitrendipine; diazepam; mianserin; ramipril; haloperidol; tetracaine; alprenolol; aceprometazine; glibenclamide; chlorophenacinone; doxepin; nimodipine; diphenhydra-

mine; cyclizine; histapyrrodine; phenylbutazone; demexiptiline; clozapine; proguanil; tri-fluperidol; medazepam; cyamemazine; bumadizone; suriclone; propranolol; acepromazine; dothiepin; dextromoramide; fenoprofen; dextropropoxyphene; loxapine; betaxolol; propafenone; promethazine; thioproperazine; methadone; amoxapine; quinupramine; opipramol; cyproheptadine; brompheniramine; mefenidramine; protriptyline; flurbiprofen; tetrazepam; zorubicin; prazepam; alimemazine; loperamide; imipramine; desipramine; levomepromazine; hydroxyzine; niflumic acid; penbutolol; fluvoxamine; pimozide; daunorubicin; indomethacin; maprotiline; tropatenine; etodolac; fluoxetine; amitriptyline; nortriptyline; tiocloamarol; diclofenac; mefloquine; trimipramine; chlorambucil; lidoflazine; ibuprofen; floctafenine; alpidem; loratadine; chlorpromazine; clomipramine; carpipramine; thioridazine; fentiazac; clemastine; mefenamic acid; fluphenazine; prochlorperazine; penfluridol; bepridil; terfenadine; trifluoperazine

REFERENCE

Tracqui,A.; Kintz,P.; Mangin,P. Systematic toxicological analysis using HPLC/DAD, *J.Forensic Sci.*, **1995**, *40*, 254–262.

SAMPLE

Matrix: blood, urine

Sample preparation: Add 1 mL whole blood or urine to Toxi-Tube A (Toxi-Lab, Irvine CA), add 3 mL water, mix by gentle inversion for 5 min, centrifuge at 1500 g for 5 min. Remove the organic layer and evaporate it to dryness under a stream of nitrogen at 40°, reconstitute the residue with 50 µL MeCN:water 50:50, vortex for 10 s, centrifuge at 7500 g for 2 min, inject a 10 (urine) or 30 (blood) µL aliquot. (The detector wavelength shown is the wavelength of maximum absorbance. This will not necessarily be the optimal wavelength for the separation. Multiple wavelengths from 200-350 nm can be scanned using a diode-array detector. Otherwise, 220 nm may be a reasonable choice for initial work. Matrix may interfere.)

HPLC VARIABLES

Guard column: 20 mm long Symmetry C18

Column: 250 × 4.6 5 µm Symmetry C8 (Waters)

Mobile phase: Gradient. A was 50 mM pH 3.8 sodium phosphate buffer. B was MeCN. A: B 85:15 for 6.5 min, 65:35 for 18.5 min, 20:80 for 3 min (step gradient), re-equilibrate at initial conditions for 7 min.

Column temperature: 30

Flow rate: 1 for 6.5 min, to 1.5 over 18.5 min, maintain at 1.5 for 3 min (re-equilibrate at 1.5 mL/min)

Injection volume: 10-30

Detector: UV 204

CHROMATOGRAM

Retention time: 20.052

KEY WORDS

whole blood

REFERENCE

Gaillard,Y.; Pépin,G. Use of high-performance liquid chromatography with photodiode-array UV detection for the creation of a 600-compound library. Application to forensic toxicology, *J.Chromatogr.A*, **1997**, *763*, 149–163.


Acepromazine

Molecular formula: C₁₉H₂₂N₂O₂S

Molecular weight: 326.46

CAS Registry No.: 61-00-7, 3598-37-6 (maleate)

Merck Index: 32

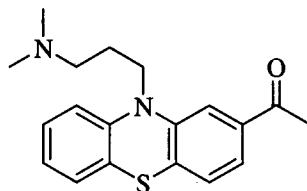

SAMPLE

Matrix: blood

Sample preparation: 2 mL Whole blood or plasma + 2 mL buffer + 5 mL chloroform: isopropanol:n-heptane 60:14:26, shake gently horizontally for 10 min, centrifuge at 2800 g for 10 min. Remove the lower organic layer and evaporate it to dryness under vacuum at 45°, reconstitute the residue in 100 µL mobile phase, centrifuge at 2800 g for 5 min, inject a 50 µL aliquot of the supernatant. (Buffer was saturated ammonium chloride solution 25% diluted with water, adjusted to pH 9.5 with 25% ammonia solution.)

HPLC VARIABLES

Column: 300 × 3.9 µm NovaPack C18

Mobile phase: MeOH:THF:buffer 65:5:30 (Buffer was 0.68 g/L (10 mM (sic)) KH₂PO₄ adjusted to pH 2.6 with concentrated orthophosphoric acid.) (At the end of each session wash the column with water for 1 h and MeOH for 1 h, re-equilibrate for 30 min.)

Column temperature: 30

Flow rate: 0.8

Injection volume: 50

Detector: UV 243

CHROMATOGRAM

Retention time: 6.77

Limit of detection: <120 ng/mL

KEY WORDS

whole blood; plasma; interferences may occur—compounds(all of which are extracted) elute in this order tenoxicam; iproniazid; methocarbamol; methotrexate; caffeine; nialamide; colchicine; cytarabine; benzoylecgonine; acetaminophen; diazoxide; dacarbazine; sulfinpyrazole; flumazenil; sulpride; morphine; atenolol; toloxatone; terbutaline; albuterol; phenobarbital; ranitidine; tiapride; phenol; chlormezanone; aspirin; metformin; ritodrine; codeine; sultopride; amisulpride; naltrexone; lisinopril; benzocaine; nizatidine; nalorphine; mephenesin; naloxone; sotalol; carteolol; procainamide; carbamazepine; bromazepam; nalbuphine; nadolol; procarbazine; dihydralazine; omeprazole; strychnine; acebutolol; glutethimide; chlorpropamide; glipizide; triazolam; prazosin; flunitrazepam; clonazepam; metoclopramide; melphalan; estazolam; tolbutamide; ephedrine; clonidine; pindolol; clobazam; minoxidil; disopyramide; nitrazepam; dextromethorphan; tofisopam; zopiclone; debrisoquine; sulindac; alprazolam; cycloguanil; lorazepam; methaqualone; ketamine; piroxicam; metoprolol; nifedipine; quinine; mephentermine; prilocaine; pentazocine; oxazepam; tiaprofenic acid; quinidine; celiprolol; ajmaline; yohimbine; lidocaine; secobarbital; viloxazine; mepivacaine; meperidine; doxylamine; labetalol; temazepam; amodiaquine; benperidol; droperidol; hydroxychloroquine; zolpidem; ketoprofen; alminoprofen; cicletanine; moclobemide; chloroquine; cocaine; timolol; nomifensine; ticlopidine; acenocoumarol; videsine; mexiletine; dipyridamole; trazodone; pipamperone; pyrimethamine; benazepril; vincristine; metapramine; chlordiazepoxide; oxprenolol; warfarin; clorazepate; flecainide; phenacyclidine; thiopental; fenfluramine; metipranolol; triprolidine; naproxen; buprenorphine; verapamil; buspirone; tianeptine; midazolam; bupivacaine; carbinoxamine; loperazolam; cetirizine; chlorpheniramine; moperone; cibenzoline; medifoxamine; astemizole; vinblastine; nicardipine; bisoprolol; diltiazem; glibornuride; reserpine; aconitine; nitrendipine; diazepam; mianserin; ramipril; haloperidol; tetracaine; alprenolol; aceprometazine; glibenclamide; chlorophenacinone; doxepin; nimodipine; diphenhydra-

mine; cyclizine; histapyrrodine; phenylbutazone; demexiptiline; clozapine; proguanil; tri-fluperidol; medazepam; cyamemazine; bumadizone; suriclone; propranolol; acepromazine; dothiepin; dextromoramide; fenoprofen; dextropropoxyphene; loxapine; betaxolol; propafenone; promethazine; thioproperazine; methadone; amoxapine; quinupramine; opipramol; cyproheptadine; brompheniramine; mefenidramine; protriptyline; flurbiprofen; tetrazepam; zorubicin; prazepam; alimemazine; loperamide; imipramine; desipramine; levomepromazine; hydroxyzine; niflumic acid; penbutolol; fluvoxamine; pimozide; daunorubicin; indomethacin; maprotiline; tropatenine; etodolac; fluoxetine; amitriptyline; nor-triptyline; tiocloamarol; diclofenac; mefloquine; trimipramine; chlorambucil; lidoflazine; ibuprofen; floctafenine; alpidem; loratadine; chlorpromazine; clomipramine; carpipramine; thioridazine; fentiazac; clemastine; mefenamic acid; fluphenazine; prochlorperazine; pen-fluridol; bepridil; terfenadine; trifluoperazine

REFERENCE

Tracqui,A.; Kintz,P.; Mangin,P. Systematic toxicological analysis using HPLC/DAD, *J.Forensic Sci.*, **1995**, *40*, 254–262.

SAMPLE

Matrix: blood, urine

Sample preparation: Add 1 mL whole blood or urine to Toxi-Tube A (Toxi-Lab, Irvine CA), add 3 mL water, mix by gentle inversion for 5 min, centrifuge at 1500 g for 5 min. Remove the organic layer and evaporate it to dryness under a stream of nitrogen at 40°, reconstitute the residue with 50 µL MeCN:water 50:50, vortex for 10 s, centrifuge at 7500 g for 2 min, inject a 10 (urine) or 30 (blood) µL aliquot. (The detector wavelength shown is the wavelength of maximum absorbance. This will not necessarily be the optimal wavelength for the separation. Multiple wavelengths from 200–350 nm can be scanned using a diode-array detector. Otherwise, 220 nm may be a reasonable choice for initial work. Matrix may interfere.)

HPLC VARIABLES

Guard column: 20 mm long Symmetry C18

Column: 250 × 4.6 5 µm Symmetry C8 (Waters)

Mobile phase: Gradient. A was 50 mM pH 3.8 sodium phosphate buffer. B was MeCN. A: B 85:15 for 6.5 min, 65:35 for 18.5 min, 20:80 for 3 min (step gradient), re-equilibrate at initial conditions for 7 min.

Column temperature: 30

Flow rate: 1 for 6.5 min, to 1.5 over 18.5 min, maintain at 1.5 for 3 min (re-equilibrate at 1.5 mL/min)

Injection volume: 10–30

Detector: UV 249.9

CHROMATOGRAM

Retention time: 10.763

KEY WORDS

whole blood

REFERENCE

Gaillard,Y.; Pépin,G. Use of high-performance liquid chromatography with photodiode-array UV detection for the creation of a 600-compound library. Application to forensic toxicology, *J.Chromatogr.A*, **1997**, *763*, 149–163.

SAMPLE

Matrix: solutions

Sample preparation: Prepare a 10 µg/mL solution in MeOH, inject a 20 µL aliquot.

HPLC VARIABLES

Column: 125 × 4.9 Spherisorb S5W silica

Mobile phase: MeOH containing 10 mM ammonium perchlorate and 1 mL/L 100 mM NaOH in MeOH, pH 6.7

Flow rate: 2

Injection volume: 20

Detector: E, LeCarbone, V25 glassy carbon electrode, + 1.2 V

CHROMATOGRAM

Retention time: 4.6

OTHER SUBSTANCES

Also analyzed: acebutolol, acetophenazine, N-acetylprocainamide, albuterol, alprenolol, amethocaine, amiodarone, amitriptyline, antazoline, atenolol, azacyclonal, bamethan, benactyzine, benperidol, benzethidine, benzocaine, benzocetamine, benzphetamine, benzquinamide, bromhexine, bromodiphenhydramine, bromperidol, brompheniramine, brompromazine, buclizine, bufotenine, bupivacaine, buprenorphine, butacaine, butethamate, chlorcyclizine, chlorpheniramine, chlorphenoxamine, chlorprenaline, chlorpromazine, chlorprothixene, cimetidine, cinchonidine, cinnarizine, clemastine, clomipramine, clonidine, cocaine, cyclazocine, cyclizine, cyclopentamine, cyproheptadine, deserpidine, desipramine, dextromoramide, dextropropoxyphene, dicyclomine, diethylcarbamazine, diethylpropion, diethylthiambutene, dihydroergotamine, dimethindene, dimethothiazine, diphenhydramine, diphenoxylate, dipipanone, diprenorphine, dipyridamole, disopyramide, dothiepin, doxapram, doxepin, doxylamine, droperidol, ephedrine, ergocornine, ergocristine, ergocristinine, ergocryptine, ergometrine, ergosine, ergosinine, ergotamine, ethopropazine, etorphine, etoxeridine, fenethazine, fenfluramine, fenoterol, fentanyl, flavoxate, fluopromazine, flupenthixol, fluphenazine, flurazepam, haloperidol, hydroxyzine, hyoscine, ibogaine, imipramine, indapamine, iprindole, isothipendyl, isoxsuprine, ketanserin, laudanosine, lidocaine, lofepramine, loxapine, maprotiline, mecamlamine, meclophenoxate, meclozine, medazepam, mephentermine, mepivacaine, meptazinol, mepyramine, mesoridazine, metaraminol, methadone, methamphetamine, methapyrilene, methdilazene, methotrimeprazine, methoxamine, methoxyphenamine, methoxypromazine, methylephedrine, methylergonovine, methysergide, metoclopramide, metopimazine, metoprolol, mianserin, morazone, nadolol, nalorphine, naloxone, naphazoline, nicotine, nifedipine, nomifensine, nortriptyline, noscapine, orphenadrine, oxeladin, oxprenolol, oxymetazolin, papaverine, pargyline, pecazine, penbutolol, pentazocine, penthienate, pericyazine, perphenazine, phenadoxone, phenampromide, phenazocine, phenbutrazate, phendimetrazine, phenelzine, phenglutarimide, phenindamine, pheniramine, phenmetrazine, phenomorphan, phenoperidine, phenothiazine, phenoxybenzamine, phentolamine, phenylephrine, phenyltoloxamine, physostigmine, piminodine, pimozide, pindolol, pipamazine, pipazethate, piperacetazine, piperidolate, pipradol, pirenzepine, piritramide, pizotifen, practolol, pramoxine, prazosin, prenylamine, prilocaine, primaquine, proadifen, procainamide, procaine, prochlorperazine, procyclidine, proheptazine, prolintane, promazine, promethazine, pronethalol, properidine, propiomazine, propranolol, prothipendyl, protriptyline, proxymetacaine, pseudoephedrine, pyrimethamine, quinidine, quinine, ranitidine, rescinnamine, sotalol, tacrine, terazosin, terbutaline, terfenadine, thenyldiamine, theophylline, thiethylperazine, thiopropazate, thiopropazine, thioridazine, thiothixene, thonzylamine, timolol, tocainide, tolpropamine, tolycaine, tranlycypromine, trazodone, trifluoperazine, trifluoperidol, trimeperidine, trimeprazine, trimethobenzamide, trimethoprim, trimipramine, tripeleminamine, triprolidine, tryptamine, verapamil, xylometazoline

REFERENCE

Jane, I.; McKinnon, A.; Flanagan, R.J. High-performance liquid chromatographic analysis of basic drugs on silica columns using non-aqueous ionic eluents. II. Application of UV, fluorescence and electrochemical oxidation detection, *J.Chromatogr.*, **1985**, 323, 191-225.

SAMPLE

Matrix: solutions

HPLC VARIABLES

Column: 250 × 4.6 Zorbax RX

Mobile phase: Gradient. A was 10 mL concentrated orthophosphoric acid and 7 mL triethylamine in 1 L water. B was 10 mL concentrated orthophosphoric acid and 7 mL triethylamine in 200 mL water, make up to 1 L with MeCN. A:B from 100:0 to 0:100 over 30 min, maintain at 0:100 for 5 min.

Column temperature: 30

Flow rate: 2

Detector: UV 210

OTHER SUBSTANCES

Also analyzed: acetaminophen, acetophenazine, albuterol, aminophylline, amitriptyline, amobarbital, amoxapine, amphetamine, amyllocaine, antipyrine, aprobarbital, aspirin, atenolol, atropine, avermectin, barbital, benzocaine, benzoic acid, benzotropine, benzphetamine, berberine, bibucaine, bromazepam, brompheniramine, buprenorphine, buspirone, butabarbital, butacaine, butethal, caffeine, carbamazepine, carbromal, chloramphenicol, chlordiazepoxide, chloroquine, chlorothiazide, chloroxylenol, chlorphenesin, chlorpheniramine, chlorpromazine, chlorpropamide, chlortetracycline, cimetidine, cinchonidine, cinchonine, clenbuterol, clonazepam, clonixin, clorazepate, cocaine, codeine, colchicine, cortisone, coumarin, cyclazocine, cyclobenzaprine, cyclothiazide, cyheptamide, cymarin, danazol, danthron, dapsone, debrisoquine, desipramine, dexamethasone, dextromethorphan, dextropropoxyphene, diamorphine, diazepam, diclofenac, diethylpropion, diethylstilbestrol, diflunisal, digitoxin, digoxin, diltiazem, diphenhydramine, diphenoxylate, diprenorphine, dipyrone, disulfiram, dopamine, doxapram, doxepin, dronabinol, ephedrine, epinephrine, epinine, estradiol, estriol, estrone, ethacrynic acid, ethosuximide, etonitazene, etorphine, eugenol, famotidine, fenbendazole, fencamfamine, fenoprofen, fenproporex, fentanyl, flubendazole, flufenamic acid, flunitrazepam, 5-fluorouracil, fluoxymesterone, fluphenazine, furosemide, gentisic acid, gitoxigenin, glipizide, glunixin, glutethimide, glybenclamide, guaiaicol, halazepam, haloperidol, hydrochlorothiazide, hydrocodone, hydrocortisone, hydromorphone, hydroxyquinoline, ibogaine, ibuprofen, iminostilbene, imipramine, indomethacin, isocarboxystyryl, isocarboxazid, isoniazid, isoproterenol, isoxsuprine, ivermectin, ketamine, ketoprofen, kynurenic acid, levorphanol, lidocaine, lorazepam, lormetazepam, loxapine, mazindol, mebendazole, meclizine, meclofenamic acid, medazepam, mefenamic acid, megestrol, mepacrine, meperidine, mephentermine, mephenytoin, mephesin, mephobarbital, mepivacaine, mescaline, mesoridazine, methadone, methamphetamine, methapyrilene, methaqualone, methazolamide, methocarbamol, methoxamine, methsuximide, methyl salicylate, methyl dopa, methyl dopamine, methylphenidate, methylprednisolone, methyltestosterone, methypylon, metoprolol, mibolerone, morphine, nadolol, nalorphine, naloxone, naltrexone, naphazoline, naproxen, nefopam, niacinamide, nicotine, niacin, nifedipine, niflumic acid, nitrazepam, norepinephrine, nortriptyline, noscapine, nylidrin, oxazepam, oxycodone, oxymorphone, oxyphenbutazone, oxytetracycline, papaverine, pargyline, pemoline, pentazocine, pentobarbital, persantine, phenacetin, phenazocine, phenazopyridine, phencyclidine, phendimetrazine, phenelzine, pheniramine, phenobarbital, phenothiazine, phensuximide, phentermine, phenylbutazone, phenylephrine, phenylpropanolamine, piperocaine, prazepam, prednisolone, primidone, probenecid, progesterone, propiomazine, propranolol, propylparaben, pseudoephedrine, puromycin, pyrilamine, pyrithione, quazepam, quinaldic acid, quinidine, quinine, ranitidine, recinnamine, reserpine, resorcinol, saccharin, albuterol, salicylamide, salicylic acid, scopolamine, scopoletin, secobarbital, strychnine, sulfacetamide, sulfadiazine, sulfadimethoxine, sulfaethidole, sulfamerazine, sulfamethazine, sulfamethoxazole, sulfanilamide, sulfapyridine, sulfasoxazole, sulindac, tamoxifen, temazepam, testosterone, tetracaine, tetracycline, tetramisole, thebaine, theobromine, theophylline, thiabendazole, thiamine, thiamylal, thiobarbituric acid, thioridazine, thiosalicylic acid, thiothixene, thymol, tolazamide, tolazoline, tobutamide, tolmetin, tranlycypromine, triamcinolone, tribenzylamine, trichloromethiazide, trifluoperazine, trihexyphenidyl, trimethoprim, tripeleminamine, triprolidine, tropacocaine, tyramine, verapamil, vincamine, warfarin, yohimbine, zoxazolamine

REFERENCE

Hill,D.W.; Kind,A.J. Reversed-phase solvent gradient HPLC retention indexes of drugs, *J.Anal.Toxicol.*, 1994, 18, 233-242.

SAMPLE

Matrix: solutions

HPLC VARIABLES

Column: 250 × 4.6 5 µm Supelcosil LC-DP (A) or 250 × 4.5 µm LiChrospher 100 RP-8 (B)
Mobile phase: MeCN:0.025% phosphoric acid:buffer 25:10:5 (A) or 60:25:15 (B) (Buffer was 9 mL concentrated phosphoric acid and 10 mL triethylamine in 900 mL water, adjust pH to 3.4 with dilute phosphoric acid, make up to 1 L.)

Flow rate: 0.6

Injection volume: 25

Detector: UV 229

CHROMATOGRAM

Retention time: 12.81 (A), 5.81 (B)

OTHER SUBSTANCES

Also analyzed: acebutolol, acetaminophen, acetazolamide, acetophenazine, albuterol, alprazolam, amitriptyline, amobarbital, amoxapine, antipyrine, atenolol, atropine, azatadine, baclofen, benzocaine, bromocriptine, brompheniramine, brotizolam, bupivacaine, buspirone, butabarbital, butalbital, caffeine, carbamazepine, cetirizine, chlorcyclizine, chlordiazepoxide, chlormezanone, chloroquine, chlorpheniramine, chlorpromazine, chlorpropamide, chlorprothixene, chlorthalidone, chlorzoxazone, cimetidine, cisapride, clomipramine, clonazepam, clonidine, clozapine, cocaine, codeine, colchicine, cyclizine, cyclobenzaprine, dantrolene, desipramine, diazepam, diclofenac, diflunisal, diltiazem, diphenhydramine, diphenidol, diphenoxylate, dipyrindamole, disopyramide, dobutamine, doxapram, doxepin, droperidol, encainide, ethidium bromide, ethopropazine, fenoprofen, fentanyl, flavoxate, fluoxetine, fluphenazine, flurazepam, flurbiprofen, fluvoxamine, furosemide, glutethimide, glyburide, guaifenesin, haloperidol, homatropine, hydralazine, hydrochlorothiazide, hydrocodone, hydromorphone, hydroxychloroquine, hydroxyzine, ibuprofen, imipramine, indomethacin, ketoconazole, ketoprofen, ketorolac, labetalol, levorphanol, lidocaine, loratadine, lorazepam, lovastatin, loxapine, mazindol, mefenamic acid, meperidine, mephénytoin, mepivacaine, mesoridazine, metaproterenol, methadone, methdilazine, methocarbamol, methotrexate, methotrimeprazine, methoxamine, methyl-dopa, methylphenidate, metoclopramide, metolazone, metoprolol, metronidazole, midazolam, moclobemide, morphine, nadolol, nalbuphine, naloxone, naphazoline, naproxen, nifedipine, nizatidine, norepinephrine, nortriptyline, oxazepam, oxycodone, oxymetazoline, paroxetine, pemoline, pentazocine, pentobarbital, pentoxifylline, perphenazine, pheniramine, phenobarbital, phenol, phenolphthalein, phentolamine, phenylbutazone, phenyltoloxamine, phenytol, pimozide, pindolol, piroxicam, pramoxine, prazepam, prazosin, probenecid, procainamide, procaine, prochlorperazine, procyclidine, promazine, promethazine, propafenone, propantheline, propiomazine, propofol, propranolol, protriptyline, quazepam, quinidine, quinine, racemethorphan, ranitidine, remoxipride, risperidone, salicylic acid, scopalamine, secobarbital, sertraline, sotalol, spironolactone, sulfinpyrazone, sulindac, temazepam, terbutaline, terfenadine, tetracaine, theophylline, thiethylperazine, thiopental, thioridazine, thiothixene, timolol, tocinide, tolbutamide, tolmetin, trazodone, triamterene, triazolam, trifluoperazine, triflupromazine, trimeprazine, trimethoprim, trimipramine, verapamil, warfarin, xylometazoline, yohimbine, zopiclone

KEY WORDS

also details of plasma extraction

REFERENCE

Koves,E.M. Use of high-performance liquid chromatography-diode array detection in forensic toxicology, *J.Chromatogr.A*, 1995, 692, 103-119.

SAMPLE**Matrix:** tissue**Sample preparation:** Condition a Sep-Pak C18 SPE cartridge with 5 mL MeOH and 5 mL water. Homogenize kidney with a kitchen grinder. Weigh out a 5 g sample and add 20 mL MeCN with continuous gentle mixing, mix vigorously on a vibromixer at 1500 rpm for 30 s, sonicate for 2 min, centrifuge at 4000 g for 5 min. Mix 7.5 mL sample extract and 40 mL 10% NaCl and add to SPE cartridge, wash with 1 mL 10 mM sulfuric acid, wash with 2 mL air, elute with 2 mL acidic MeCN. Place eluate in a washed tube and evaporate to 300 μ L at 70° under a stream of nitrogen, mix gently, add 1 mL n-hexane, mix on a vibromixer for 30 s, centrifuge at 2000 g, inject a 50 μ L aliquot of the aqueous phase. (Acidic MeCN was 1 mL 50 mM sulfuric acid and 100 mL MeCN. The washed tube was prepared by rinsing with concentrated ammonia, water, and acetone and drying under a stream of nitrogen.)

HPLC VARIABLES**Guard column:** 10 \times 2.1 37-50 μ m Bondapak C18**Column:** 300 \times 3.9 Bondapak C18**Mobile phase:** MeCN:water 55:45 containing 2.46 g/L anhydrous sodium acetate, pH adjusted to 6.5 with acetic acid**Flow rate:** 1.2**Injection volume:** 50**Detector:** UV 240

CHROMATOGRAM**Retention time:** 12**Limit of detection:** 2 ng/g

OTHER SUBSTANCES**Extracted:** azaperol, carazolol, xylazine, azaperone, haloperidol, propiomazine, chlorpromazine

KEY WORDS

SPE; pig; kidney

REFERENCE

Keukens,H.J.; Aerts,M.M.L. Determination of residues of carazolol and a number of tranquilizers in swine kidney by high-performance liquid chromatography with ultraviolet and fluorescence detection, *J.Chromatogr.*, **1989**, 464, 149-161.

SAMPLE**Matrix:** tissue**Sample preparation:** Condition a Bond-Elut C18 SPE cartridge with 5 mL MeOH and 5 mL water. Cut pig kidney or liver into small pieces and homogenize. 5 g Homogenate + 10 mL MeCN, shake, vortex for 30 s, sonicate for 3 min, vortex for 30 s, sonicate for 3 min, centrifuge at 10000 g for 20 min. Add 7.5 mL supernatant + 40 mL 10% NaCl to the SPE cartridge at about 1 mL/min, do not allow cartridge to dry out, wash with 850 μ L 10 mM sulfuric acid, dry with air, elute with 3.5 mL acidic MeCN. Evaporate the eluate to dryness under a stream of nitrogen at 50°, reconstitute the residue in 300 μ L 10 mM sulfuric acid, vortex briefly, add 1 mL hexane, vortex for 30 s, centrifuge at 2000 g for 5 min, inject an aliquot of the aqueous layer. (Acidic MeCN was 1 mL 50 mM sulfuric acid in 100 mL MeCN.)

HPLC VARIABLES**Guard column:** Hypersil 5 μ m SAS C1**Column:** 250 mm long 5 μ m Hypersil SAS C1**Mobile phase:** MeCN:water 50:50 containing 0.77 g/L ammonium acetate**Flow rate:** 2

Detector: E, ESA Model 5100A Coulochem, first electrode +0.4 V, second electrode (which was monitored) +0.7 V, Model 5020 guard cell after pump but before injector at +0.75 V

CHROMATOGRAM

Retention time: 18

Limit of detection: 2 ng/g

OTHER SUBSTANCES

Extracted: azaperol, azaperone, carazolol, xylazine, haloperidol, propriomazine, chlorpromazine

KEY WORDS

SPE; pig; kidney; liver

REFERENCE


Rose,M.D.; Shearer,G. Determination of tranquilisers and carazolol residues in animal tissue using high-performance liquid chromatography with electrochemical detection, *J.Chromatogr.*, **1992**, 624, 471-477.

Aceprometazine

Molecular formula: C₁₉H₂₂N₂O₂S

Molecular weight: 326.46

CAS Registry No.: 13461-01-3

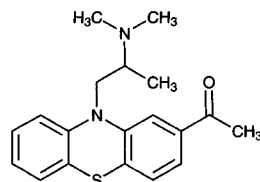

SAMPLE

Matrix: blood

Sample preparation: 2 mL Whole blood or plasma + 2 mL buffer + 5 mL chloroform: isopropanol:n-heptane 60:14:26, shake gently horizontally for 10 min, centrifuge at 2800 g for 10 min. Remove the lower organic layer and evaporate it to dryness under vacuum at 45°, reconstitute the residue in 100 µL mobile phase, centrifuge at 2800 g for 5 min, inject a 50 µL aliquot of the supernatant. (Buffer was saturated ammonium chloride solution 25% diluted with water, adjusted to pH 9.5 with 25% ammonia solution.)

HPLC VARIABLES

Column: 300 × 3.9 4 µm NovaPack C18

Mobile phase: MeOH:THF:buffer 65:5:30 (Buffer was 0.68 g/L (10 mM (sic)) KH₂PO₄ adjusted to pH 2.6 with concentrated orthophosphoric acid.) (At the end of each session wash the column with water for 1 h and MeOH for 1 h, re-equilibrate for 30 min.)

Column temperature: 30

Flow rate: 0.8

Injection volume: 50

Detector: UV 242

CHROMATOGRAM

Retention time: 6.20

Limit of detection: <120 ng/mL

KEY WORDS


whole blood; plasma; interferences may occur—compounds(all of which are extracted) elute in this order tenoxicam; iproniazid; methocarbamol; methotrexate; caffeine; nialamide; colchicine; cytarabine; benzoylecgonine; acetaminophen; diazoxide; dacarbazine; sulfinpyrazole; flumazenil; sulpride; morphine; atenolol; tolaxatone; terbutaline; albuterol; phenobarbital; ranitidine; tiapride; phenol; chlormezanone; aspirin; metformin; ritodrine; codeine; sultopride; amisulpride; naltrexone; lisinopril; benzocaine; nizatidine; nalorphine; mephenesin; naloxone; sotalol; carteolol; procainamide; carbamazepine; bromazepam; nalbuphine; nadolol; procarbazine; dihydralazine; omeprazole; strychnine; acebutolol; glutethimide; chlorpropamide; glipizide; triazolam; prazosin; flunitrazepam; clonazepam; metoclopramide; melphalan; estazolam; tolbutamide; ephedrine; clonidine; pindolol; clobazam; minoxidil; disopyramide; nitrazepam; dextromethorphan; tofisopam; zopiclone; debrisoquine; sulindac; alprazolam; cycloguanil; lorazepam; methaqualone; ketamine; piroxicam; metoprolol; nifedipine; quinine; mephentermine; prilocaine; pentazocine; oxazepam; tiaprofenic acid; quinidine; celiprolol; ajmaline; yohimbine; lidocaine; secobarbital; viloxazine; mepivacaine; meperidine; doxylamine; labetalol; temazepam; amodiaquine; benperidol; droperidol; hydroxychloroquine; zolpidem; ketoprofen; alminoprofen; cicletanine; mocllobemide; chloroquine; cocaine; timolol; nomifensine; ticlopidine; acenocoumarol; vindsine; mexiletine; dipyrindamole; trazodone; pipamperone; pyrimethamine; benazepril; vincristine; metapramine; chlordiazepoxide; oxprenolol; warfarin; clorazepate; flecainide; phencyclidine; thiopental; fenfluramine; metipranolol; triprolidine; naproxen; buprenorphine; verapamil; buspirone; tianeptine; midazolam; bupivacaine; carbinoxamine; loperazolam; cetirizine; chlorpheniramine; moperone; cibenzone; medifoxamine; astemizole; vinblastine; nicardipine; bisoprolol; diltiazem; glibornuride; reserpine; aconitine; nitrendipine; diazepam; mianserin; ramipril; haloperidol; tetracaine; alprenolol; aceprometazine; glibenclamide; chlorophenacinone; doxepin; nimodipine; diphenhydramine; cyclizine; histapyrodine; phenylbutazone; demexiptiline; clozapine; proguanil; tri-

fluperidol; medazepam; cyamemazine; bumadizone; suriclone; propranolol; acepromazine; dothiepin; dextromoramide; fenoprofen; dextropropoxyphene; loxapine; betaxolol; propafenone; promethazine; thioproperazine; methadone; amoxapine; quinupramine; opipramol; cyproheptadine; brompheniramine; mefenidramine; protriptyline; flurbiprofen; tetrazepam; zorubicin; prazepam; alimemazine; loperamide; imipramine; desipramine; levomepromazine; hydroxyzine; niflumic acid; penbutolol; fluvoxamine; pimozone; daunorubicin; indomethacin; maprotiline; tropatenine; etodolac; fluoxetine; amitriptyline; nortriptyline; tiocloamarol; diclofenac; mefloquine; trimipramine; chlorambucil; lidoflazine; ibuprofen; floctafenine; alpidem; loratadine; chlorpromazine; clomipramine; carpipramine; thioridazine; fentiazac; clemastine; mefenamic acid; fluphenazine; prochlorperazine; penfluridol; bepridil; terfenadine; trifluoperazine

REFERENCE

Tracqui,A.; Kintz,P.; Mangin,P. Systematic toxicological analysis using HPLC/DAD, *J.Forensic Sci.*, **1995**, *40*, 254–262.

Acesulfame

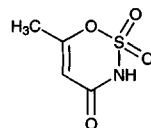

Molecular formula: C₄H₅NO₄S

Molecular weight: 163.15

CAS Registry No.: 33665-90-6

SAMPLE

Matrix: beverages, sweetener

Sample preparation: Sweetener. Dissolve 30 mg powdered tabletop sweetener in water and dilute to 25 mL, filter (0.2 µm PTFE). Beverages. Dilute fruit juice 1:10 with water. Degas carbonated beverages in a ultrasonic bath for 5 min, dilute 1:10 with water, filter. Inject a 50 µL aliquot.

HPLC VARIABLES

Guard column: 50 × 4 Dionex IonPak AG4A-SC

Column: 250 × 4 Dionex IonPak AS4A-SC

Mobile phase: Gradient. A was 1 mM sodium carbonate. B was 12.5 mM sodium carbonate. A:B 100:0 for 4.5 min, from 100:0 to 0:100 over 1 min, maintain at 0:100 for 22.5 min, from 0:100 to 100:0 over 0.1 min

Flow rate: 1

Injection volume: 50

Detector: UV 190 for 6 min, UV 206 22 min, then UV 190; Conductivity, Dionex ED40 in conductivity mode preceded by a Dionex ASRS-I suppressor (external water mode, 300 mA)

CHROMATOGRAM

Retention time: 14

Limit of detection: 44 ng/mL (UV), 230 ng/mL (conductivity)


OTHER SUBSTANCES

Simultaneous: aspartame, saccharin

REFERENCE

Chen,Q.-C.; Mou,S.-.; Liu,K.-.; Yang,Z.-.; Ni,Z.-. Separation and determination of four artificial sweeteners and citric acid by high-performance liquid chromatography, *J.Chromatogr.A*, **1997**, 771, 135-143.

Acetaminophen

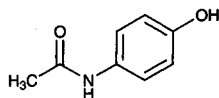

Molecular formula: $C_8H_9NO_2$

Molecular weight: 151.16

CAS Registry No.: 103-90-2

Merck Index: 45

SAMPLE

Matrix: bile, blood, gastric contents, tissue, urine

Sample preparation: Chop 5-g tissue and homogenize (Ultra Turrax T25) at 8500, 9500, 13500, 20500, and 24000 rpm for 1 min each. Add homogenate to 20 mL water. Dilute blood, urine, gastric contents, and bile four times with water. Mix 4 mL sample with 10 μ L 200 μ g/mL IS and 1 mL pH 7.4 phosphate buffer, vortex briefly, add 4 mL diethyl ether and mix for 15 min (Spiramix 10, Denley, UK). Separate the organic layer, add 4 mL diethyl ether to extraction sample, mix. Evaporate combined organic layers to dryness under a stream of dry air at 50°. Purify extracts by partition between 1 mL MeCN and 2 mL heptane, separate MeCN layer, evaporate it to dryness, reconstitute the residue in 1 mL MeOH and inject a 20 μ L aliquot of the solution.

HPLC VARIABLES

Guard column: 20 \times 4.6 5 μ m Apex II ODS

Column: 150 \times 4.6 5 μ m Apex II ODS

Mobile phase: MeCN:acetic acid:water 10:5:85

Flow rate: 1

Injection volume: 20

Detector: UV 255

CHROMATOGRAM

Retention time: 3.3

Internal standard: 2-acetoamidophenol (6.4)

Limit of quantitation: 2 μ g/mL

KEY WORDS

liver; lung; muscle; urine; pericardial fluid

REFERENCE

Pounder,D.J.; Adams,E.; Fuke,C.; Langford,A.M. Site to site variability of postmortem drug concentrations in liver and lung, *J.Forensic Sci.*, **1996**, *41*, 927-932.

SAMPLE

Matrix: blood

Sample preparation: 200 μ L Plasma + 200 μ L 1.0 M perchloric acid, centrifuge. Add 200 μ L 700 mM potassium phosphate, cool on ice for 30 min, centrifuge for 5 min. Inject an aliquot of the supernatant.

HPLC VARIABLES

Column: 100 \times 4.6 C8 Rainin "Short One"

Mobile phase: MeCN:buffer 4:96 (Buffer was 25 mM phosphate with 0.5% acetic acid, pH 3.1)

Flow rate: 1.0

Detector: E, HP 1049A amperometric detector, 700 mV, oxidation mode

CHROMATOGRAM

Retention time: 6.5

KEY WORDS

plasma

REFERENCE

Sarich,T.; Kalhorn,T.; Magee,S.; Al-sayegh,F.; Adams,S.; Slattery,J.; Goldstein,J.; Nelson,S.; Wright,J.
The effect of omeprazole pretreatment on acetaminophen metabolism in rapid and slow metabolizers
of S-mephenytoin, *Clin.Pharmacol.Ther.*, **1997**, 62, 21–28.

SAMPLE

Matrix: blood

Sample preparation: Inject a 5-20 μL aliquot of serum directly.

HPLC VARIABLES

Column: 100 \times 4.6 5-10 μm Silicalite (by sieving Silicalite, 3M Co.(?))

Mobile phase: MeCN:20 mM pH 6.9 phosphate buffer 5:95 (A) or Gradient. MeCN:20 mM
pH 6.9 phosphate buffer from 5:95 to 20:80 over 2 min, to 25:75 over 2 min, to 30:70 over
4 min, to 50:50 over 2 min, maintain at 50:50 for 10 min (B)

Flow rate: 1

Injection volume: 5 (A), 20 (B)

Detector: UV 254

CHROMATOGRAM

Retention time: 4.99 (A), 6 (B)

OTHER SUBSTANCES

Extracted: barbital (B), carbamazepine (B), phenobarbital (B), phenytoin (B), primidone
(B), sulfapyridine (B), theophylline (A),

Also analyzed: metabolites

KEY WORDS

serum

REFERENCE

Ambrose,D.L.; Fntz,J.S. High-performance liquid chromatographic determination of drugs and metab-
olites in human serum and urine using direct injection and a unique molecular sieve,
J.Chromatogr.B, **1998**, 709, 89–96.

SAMPLE

Matrix: blood, urine

Sample preparation: Add 1 mL whole blood or urine to Toxi-Tube A (Toxi-Lab, Irvine CA),
add 3 mL water, mix by gentle inversion for 5 min, centrifuge at 1500 g for 5 min. Remove
the organic layer and evaporate it to dryness under a stream of nitrogen at 40°, recon-
stitute the residue with 50 μL MeCN:water 50:50, vortex for 10 s, centrifuge at 7500 g
for 2 min, inject a 10 (urine) or 30 (blood) μL aliquot. (The detector wavelength shown is
the wavelength of maximum absorbance. This will not necessarily be the optimal wave-
length for the separation. Multiple wavelengths from 200-350 nm can be scanned using
a diode-array detector. Otherwise, 220 nm may be a reasonable choice for initial work.
Matrix may interfere.)

HPLC VARIABLES

Guard column: 20 mm long Symmetry C18

Column: 250 \times 4.6 5 μm Symmetry C8 (Waters)

Mobile phase: Gradient. A was 50 mM pH 3.8 sodium phosphate buffer. B was MeCN. A:
B 85:15 for 6.5 min, 65:35 for 18.5 min, 20:80 for 3 min (step gradient), re-equilibrate at
initial conditions for 7 min.

Column temperature: 30

Flow rate: 1 for 6.5 min, to 1.5 over 18.5 min, maintain at 1.5 for 3 min (re-equilibrate at 1.5 mL/min)

Injection volume: 10-30

Detector: UV 200.5

CHROMATOGRAM

Retention time: 5.592

KEY WORDS

whole blood

REFERENCE

Gaillard, Y.; Pépin, G. Use of high-performance liquid chromatography with photodiode-array UV detection for the creation of a 600-compound library. Application to forensic toxicology, *J.Chromatogr.A*, **1997**, 763, 149-163.

SAMPLE

Matrix: formulations

Sample preparation: Add 50 mL of mobile phase to 0.5 g of sample and swirl to aid dissolution. Dilute to 100 mL with mobile phase. Dilute 1:10, filter (0.22 µm nylon). Inject a 10 µL aliquot.

HPLC VARIABLES

Column: 100 × 2.1 5 µm Hypersil ODS

Mobile phase: MeCN:water:triethylamine:acetic acid 5.5:94.1:0.2:0.2

Flow rate: 1.5

Injection volume: 10

Detector: UV 254

CHROMATOGRAM

Retention time: 0.7

OTHER SUBSTANCES

Simultaneous: aspirin, caffeine

KEY WORDS

powder

REFERENCE

Ferguson, G.K. Quantitative HPLC analysis of an analgesic/caffeine formulation: Determination of caffeine, *J.Chem.Educ.*, **1998**, 75, 467-469.

SAMPLE

Matrix: formulations

Sample preparation: Weigh 500 mg homogenized analgesic powder, transfer to 100 mL volumetric flask, add ca. 50 mL mobile phase, swirl and dilute to volume with mobile phase. Dilute an aliquot of this solution 1:10 with mobile phase, filter (0.20 µm Cameo nylon filter, MSI, Westboro, MA) an aliquot, inject an aliquot of the filtrate.

HPLC VARIABLES

Column: 100 × 2.1 5 µm Hypersil ODS

Mobile phase: MeCN:triethylamine:acetic acid:water 5.5:0.2:0.2:94.1 (Prepare mobile phase as follows. Mix 110 mL MeCN, 4 mL triethylamine, 4 mL glacial acetic acid and make up to 2 L with water.)

Flow rate: 1.5

Injection volume: 10

Detector: UV 254

CHROMATOGRAM

Retention time: 0.7

OTHER SUBSTANCES

Extracted: aspirin, caffeine

Noninterfering: salicylic acid

KEY WORDS

powder

REFERENCE

Ferguson, G.K. Quantitative HPLC analysis of an analgesic/caffeine formulation: Determination of caffeine, *J.Chem.Educ.*, **1998**, 75, 467–469.

SAMPLE

Matrix: formulations

Sample preparation: Weigh out powdered sample containing 51 mg acetaminophen, add 80 mL MeOH, sonicate for 10 min, dilute to 100 mL with MeOH, centrifuge. Remove a 5 mL aliquot of the supernatant and add it to 1 mL 2 mg/mL resorcinol, add 2 mL MeOH, make up to 20 mL with 50 mM pH 3.0 triethylamine phosphate, inject an aliquot.

HPLC VARIABLES

Column: 150 × 3.2 5 µm Hypersil ODS

Mobile phase: THF:50 mM pH 3.0 triethylamine phosphate 12:88

Flow rate: 0.6

Injection volume: 20

Detector: UV 275

CHROMATOGRAM

Retention time: 15

Internal standard: resorcinol (9)

OTHER SUBSTANCES

Simultaneous: aspirin (post-column irradiation gives an increase in peak height), caffeine, propyphenazone

REFERENCE

Di Pietra, A.M.; Gatti, R.; Andrisano, V.; Cavrini, V. Application of high-performance liquid chromatography with diode-array detection and on-line post-column photochemical derivatization to the determination of analgesics, *J.Chromatogr.A*, **1996**, 729, 355–361.

SAMPLE

Matrix: solutions

Sample preparation: Inject a 20 µL aliquot.

HPLC VARIABLES

Column: 250 × 4.6 Spherisorb S10-ODS2

Mobile phase: MeOH:20 mM pH 4.0 acetate buffer 13:87

Flow rate: 1.3

Detector: UV 243, UV 254

CHROMATOGRAM

Internal standard: 2-acetamidophenol

REFERENCE

Galia, E.; Nicolaides, E.; Hörter, D.; Löbenberg, R.; Reppas, C.; Dressman, J.B. Evaluation of various dissolution media for predicting in vivo performance of class I and II drugs, *Pharm.Res.*, **1998**, 15, 698–705.

SAMPLE**Matrix:** urine**Sample preparation:** Unhydrolyzed urine. 200 μ L Urine + 600 μ L 200 mM pH 5.0 sodium acetate/acetic acid buffer, filter through a 10 kDa molecular weight cut-off membrane (Alltech) while centrifuging. Inject a 5 μ L aliquot of the filtrate. Hydrolyzed urine. 500 μ L Urine + 500 μ L 125 mM pH 5.0 sodium acetate buffer containing 40 μ L β -glucuronidase/sulfatase, incubate overnight at 37°. Filter a 500 μ L aliquot through a 10 kDa molecular weight cut-off membrane (Alltech) while centrifuging. Inject a 10 μ L aliquot of the filtrate.

HPLC VARIABLES**Column:** 100 mm long 3 μ m Microsorb-MV C8**Mobile phase:** Gradient. MeCN:MeOH:buffer from 0:5:95 to 0:9:91 over 10 min, to 0:20:80 over 10 min, to 5:20:75 over 1 min, to 30:20:50 over 1 min. (Buffer was 25 mM pH 3.4 dibasic ammonium phosphate-acetate, 0.15% (sic).)**Injection volume:** 5-10**Detector:** UV 254

CHROMATOGRAM**Retention time:** 16-17

OTHER SUBSTANCES**Extracted:** metabolites

REFERENCE

Sarich,T.; Kalhorn,T.; Magee,S.; Al-sayegh,F.; Adams,S.; Slattery,J.; Goldstein,J.; Nelson,S.; Wright,J.
The effect of omeprazole pretreatment on acetaminophen metabolism in rapid and slow metabolizers of S-mephenytoin, *Clin.Pharmacol.Ther.*, **1997**, 62, 21-28.